


PRESSURIZED SOLUTION DISTRIBUTION SYSTEM (PSDS)

Operation and Maintenance Manual

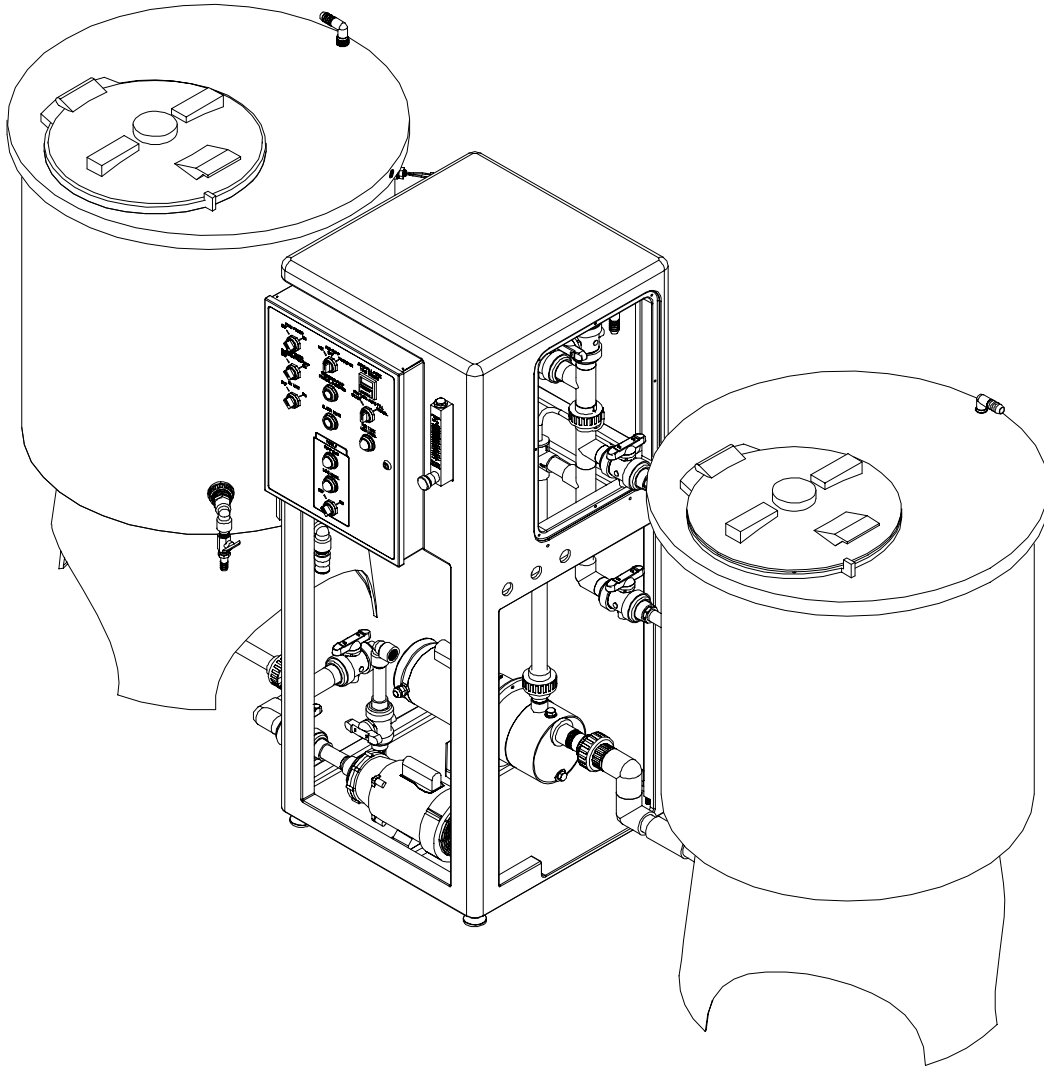

**PRESSURIZED SOLUTION DISTRIBUTION SYSTEM
(PSDS)
Operation and Maintenance Manual**

TABLE OF CONTENTS

			Page
CHAPTER ONE:		GENERAL INFORMATION	1-1
Section	1.1	The Manual	1-3
	1.2	Safety Summary	1-3
	1.3	Applications	1-4
	1.4	Contraindications	1-4
	1.5	Environmental Considerations	1-5
	1.6	Theory of Operation	1-6
	1.7	Fluid Components	1-6
	1.8	Monitors and Controls	1-8
	1.9	Alarms	1-10
	1.10	Options	1-11
	1.11	Specifications	1-13
	1.12	Disposal	1-13
	1.13	Electromagnetic Interference	1-14
	1.14	Service Assistance	1-14
	1.15	Return Materials Authorization (RMA) Procedure	1-15
	1.16	Symbols and Abbreviations	1-16
CHAPTER TWO:		INSTALLATION GUIDELINES AND INSTRUCTIONS	2-1
Section	2.1	Installation Guidelines	2-3
	2.2	PSDS Installation Requirements	2-3
	2.3	Distribution Loop Installation Requirements	2-4
	2.4	Patient Station Instructions	2-4
	2.5	Installation Instructions	2-5
	2.6	Bicarb Distribution Loop Connections	2-6
	2.7	Acid Pump Controller Connections	2-7
	2.8	Acid Distribution Loop Connections	2-9
	2.9	PSDS Status Monitor Installation	2-9

CHAPTER THREE:	SYSTEM OPERATION	3-1	
Section	3.1	Daily Start-Up	3-3
	3.2	Acid Distribution Start-Up Procedure	3-3
	3.3	Mix Tank Fill Procedure	3-4
	3.4	Bicarb Solution Mix Procedure	3-5
	3.5	Bicarb Solution Transfer Procedure	3-6
	3.6	Bicarb Solution Distribution Procedure	3-7
	3.7	Bicarb Distribution Tank & Loop Drain Procedure	3-8
	3.8	PSDS Rinse Procedure	3-12
	3.9	PSDS Disinfection/De-Scaling Procedure	3-15
	3.10	PSDS Status Monitor	3-21
CHAPTER FOUR:	ROUTINE MAINTENANCE	4-1	
Section	4.1	Routine Maintenance of Solution Distribution System	4-3
	4.2	Specific Maintenance Procedures and Instructions	4-6
	4.2.1	Bicarb System Disinfection or De-scalification Procedure	4-6
	4.2.2	Recommended Bacterial Monitoring Procedure	4-6
	4.2.3	Leak Repair	4-7
	4.2.4	Bicarb Mix Pump Maintenance	4-8
	4.2.5	PSDS Status Monitor Maintenance	4-10

CHAPTER FIVE:	TROUBLESHOOTING	5-1
	Table of Contents	5-3
Section	5.1 Alarm System	5-5
	5.2 Auto-Fill Mode Time Pre-set	5-11
	5.3 The Bicarb Mixing System	5-11
	5.4 Distribution Loop Problems (Acid/Bicarb)	5-13
	5.5 Leaks	5-14
	5.6 Rinse Pump	5-15
	5.7 Main Power Switch	5-15
	5.8 Programmable Logic Controller (PLC)	5-16
	5.9 Inlet Solenoid Valve	5-17

CHAPTER SIX: DRAWINGS 6-1

CHAPTER SEVEN: SPARE PARTS LIST 7-1

APPENDIX A: TECHNOTES/LOG SHEETS

Mar Cor is a registered trademark of Evoqua Water Technologies LLC, in some countries.

MINNCARE HD® is a register trademark of Medivators

©2024 Evoqua Water Technologies LLC. All rights reserved.

This page intentionally left blank.

**Pressurized Solution Distribution System
PSDS**

**CHAPTER ONE:
GENERAL INFORMATION**

This page intentionally left blank.

1.1 The Manual

This manual has been prepared to provide the operator with information and instructions regarding the installation, use, maintenance, and troubleshooting of the Pressurized Solution Distribution System (PSDS).

CAUTION:

When used as a medical device, Federal law restricts this device to sale by or on the order of a physician. Per CFR 801.109 (b)(1).

The manual has been written in narrative form supplemented with schematics and drawings for clarification. The operator can perform most procedures mentioned in this manual. Any exceptions will be clearly identified by a qualifying statement.

1.2 Safety Summary

Words in **BOLD CAPITAL** letters are used to identify key safety or qualifying statements. A list of all symbols and abbreviations is located at the end of this chapter.

This safety summary does not contain all of the safety statements in the manual. Other safety statements are included within the manual text and are enhanced and defined as follows:

WARNING:

Statements identifying conditions or practices that could result in personal injury or loss of life.

CAUTION:

Statements identifying conditions or practices that could result in equipment or other property damage.

NOTE:

Statements that provide further clarification.

READ THIS MANUAL:

Prior to operating or servicing this device, this manual must be read and understood. Keep this and other associated manuals for future reference and for new operators or qualified service personnel. A note sheet is provided at the end of each chapter for operators to make notations that may be valuable to other users.

USE PROPER POWER CONNECTIONS:

Use proper wiring and connection methods to satisfy hospital electrical codes.

DO NOT REMOVE COVERS OR PANELS:

To avoid electrical shock hazard, do not remove covers or panels when power is supplied to the device. Do not operate the device when covers or panels are removed.

SHOCK HAZARD:

Connect this device to a proper ground connection in accordance with the National Electrical Code. DO NOT under any circumstances remove the ground wire or ground prong from any power plug. DO NOT use an extension cord with this equipment.

DEVICE LABELING:

Do not, under any circumstances, remove any Caution, Warning or any other descriptive labels from the devices until the conditions warranting the label are eliminated.

TRANSPORTING:

When lifting or carrying the PSDS unit use at least 4 persons or properly rated lifting equipment. Ensure all proper safety equipment is used when moving the PSDS.

DO NOT OPERATE IN A FLAMMABLE ATMOSPHERE:

To avoid fire or explosion, do not operate this device in an explosive environment or near flammable anesthetics.

1.3 Applications

The Pressurized Solution Distribution System (PSDS) is to be used in a hemodialysis facility for mixing and distribution of sodium bicarbonate solution and the distribution of acid concentrate solution to remote points of use where they are used in hemodialysis.

WARNING:

The PSDS is not equipped with solution monitoring devices for conductivity, pH, or other chemically founded parameters. The operator is responsible for following the testing recommendations of the solution, chemical, or product manufacturer.

NOTE:

The PSDS is a solution distribution system and does not perform proportioning functions.

CAUTION:

If the PSDS is used in a manner not specified by the manufacturer the protective features of the unit might be impaired.

1.4 Contraindications

The Pressurized Solution Distribution System is not designed, sold, or intended to be used outside of the device specifications and limitations, as outlined in this manual and other related materials.

1.5 Environmental Considerations

Prior to the installation of your PSDS, it will be necessary to provide utilities and create an environment suitable for the trouble free operation. Performance will be affected by compliance.

POWER:

The PSDS operates on 115V single-phase power. Histories of power failure, power surges, and low line voltages should be noted and reported to the manufacturer or their agent as they may create adverse conditions for the operation of equipment. A backup power source or line conditioner may be required for uninterrupted operation.

DRAIN:

The PSDS requires a drain outlet. The drain must have a minimum capacity of twenty gallons per minute of continuous flow. The maximum height for the drain is six inches. The PSDS drain(s) should be routed separately from other equipment to a drain.

TEMPERATURE:

The PSDS should be located in an environment that will protect it from freezing or excessive heat.

NOISE:

It is advisable to locate the system away from the patient area. Hard walls will reflect noise and may make the unit seem louder.

WATER:

The PSDS requires water that meets AAMI or local hemodialysis water quality standards. Additionally, it is necessary to have adequate flow rates and pressures. The minimum flow rate required for operation of the PSDS is ½ gallon per minute.

The supplied water temperature affects the bicarbonate mixing ability. Follow the bicarbonate manufacturer's recommendations for preferred water temperature.

NOTE:

Three to five gallons per minute flow rate is recommended. Lower flow rates will prolong the filling, rinsing and disinfection procedures of the mix and distribution tanks.

NOTE:

Low water temperature may cause difficulty in mixing and/or will result in extended mixing times.

1.6 Theory of Operation

The Pressurized Solution Distribution System (PSDS) provides semi-automatic mixing of a bicarbonate concentrate powder with supplied AAMI quality water and distributes this solution along with up to three acid concentrates, depending on the model, to hemodialysis patient stations.

The bicarbonate solution is thoroughly mixed in a mix tank before being transferred to a distribution tank for circulation through the distribution loop. The distribution and mix tanks, pumps, plumbing, and controls are connected or located next to the PSDS frame.

The operational pressure in the distribution loop must not be allowed to exceed the dialysis machine manufacturer's maximum pressure; therefore, the organization recommends the use of pressure regulation at each point of use to prevent these issues.

The PSDS is available in models that can control the distribution of acid concentrate from a bulk storage tank to remote points of use. The acid distribution pumps are located near the bulk tanks with controls centralized on the PSDS frame.

1.7 Fluid Components

Inlet Solenoid Valve

The inlet solenoid valve serves to control the supply of water to the mix tank and rinsing procedure.

Flow Meter

The flow meter includes an integral needle valve to provide visual flow rate indication and adjustment.

Mix Tank

The mix tank is a high purity polyethylene tank with a conical shaped bottom and molded volume markings. The mix tank includes the following components:

- The mix tank is equipped with a spray head for rinsing and disinfecting internal surfaces.
- A sampling port is located at the bottom of the mix tank for testing/sampling purposes.
- A hinged cover is provided on top of the mix tank for easy access.
- An eductor (mix nozzle) is located at the bottom of the mix tank to aid in mixing.
- A carbon vent filter connected to the top of the mix (and distribution) tanks inhibits dust or other debris from entering.

Drains

Two drain valves are provided to allow discharge of unused solution or rinse water. There is one drain on each of the mix and distribution tanks.

Mix Pump

The mix pump provides for mixing of the bicarbonate powder into solution and transfer of the mixed solution from the mix tank to the distribution tank. The mix pump is controlled with a timer that will shutdown the pump after 10 minutes of operation to minimize the potential of over mixing. If more than 10 minutes of mixing is required, the hand switch will need to be turned to OFF and back to MIX to resume mixing.

Flow Switch

A flow switch is located on the output side of the mix pump. The flow switch prevents damage to the mix pump by automatically turning the mix pump off during insufficient flow situations. An audible alarm will also sound.

Distribution Tank

The distribution tank is a high purity polyethylene tank with a conical shaped bottom and molded volume markings. The distribution tank includes the following components:

- The distribution tank includes float switches to communicate the high-level and low-level conditions to the PSDS controls.
- The distribution tank is equipped with a spray head for rinsing and disinfecting internal surfaces.
- A jug fill valve is located near the bottom of the distribution tank for filling containers with mixed bicarbonate solution.
- A hinged cover is provided on top of the distribution tank for easy access.
- A carbon vent filter connected to the top of the distribution (and mix) tanks inhibits dust or other debris from entering.
- A valve is located near the return of the distribution loop to allow sampling of actual loop conditions.

Distribution Pump

The distribution pump(s) provide for delivery of the bicarbonate and acid solutions from the distribution tank or bulk acid tank(s) into and throughout the distribution loop. The organization recommends that over pressure protection is installed at each use point.

1.8 Monitors and Controls

Main Power Switch

The Main Power switch is located on the front of the control panel and must be turned to the ON position to operate the PSDS.

Mix Pump Switch

The Mix Pump switch is a three-position switch that controls power to the mix/transfer pump. With the switch in the TRANSFER position, the distribution tank high-level alarm and mix tank low flow alarms are both engaged and either can stop the pump operation.

Auto Fill Mode Timer

The Auto Fill Mode timer is used along with the Mix Tank Auto Fill switch and inlet water solenoid to fill the mix tank for a predetermined amount of time. The timer mode must be set to “E” and “m” (minutes) with user adjustable time settings. The time setting will control the duration the inlet water solenoid opens to supply the mix tank with water.

Mix Tank Auto Fill Switch

The Mix Tank Auto Fill switch is a momentary switch that initiates automatic or manual filling of the mix tank. When the switch is moved to the START position, the fill timer and solenoid are activated, filling the mix tank for the duration set by the timer. When the switch is held in the MANUAL position, the solenoid is activated, allowing water to flow into the mix tank as long as the switch remains held in this position.

Alarm Mute Switch

The Alarm Mute switch is a push button switch that when pressed will silence the audible alarm. The alarm will not sound again if the alarm condition is removed during this period with the exception of the Distribution Tank Low Level alarm that does require the resetting of the Distribution Tank Low Level hand switch.

Auto Fill Flow Meter

The flow meter scale allows visual setting between 0.5 and 5.0 GPM. The Auto Fill flow meter includes an integral needle valve to provide visual flow rate indication and adjustment.

Programmable Logic Controller (PLC)

The programmable logic controller (PLC) is an input/output logic device that does not perform any calculations. Software is loaded into the PLC during the assembly and does not require operator interface once installed. The PLC continuously senses for inputs from hand switches, float switches, and the mix pump flow switch. PLC outputs control operation of the mix pump, auto fill timer, and audible alarms.

Distribution Tank Low Level Alarm Switch

The Distribution Tank Low Level Alarm hand switch controls the function of the distribution tank low-level alarm. The alarm is enabled when the switch is in the ON position. Moving the switch to the BYPASS/RESET position disables the alarm and illuminates the switch to provide additional user awareness of the alarm bypass.

Distribution Pump Switch

The Distribution Pump switch is a two-position switch that controls power to the distribution pump. Moving the Distribution Pump switch to the ON position initiates operation of the Distribution pump and illuminates the Distribution Pump Run light.

Acid Tank Pump Switch

The Acid Tank Pump switch is a two-position switch that controls power to the acid tank pump. Moving the Acid Tank pump switch to the ON position initiates operation of the Acid Tank pump and illuminates the Acid Tank Pump Run light. Up to three acids can be operated.

Distribution Pump Run Light

The green Distribution Pump Run light illuminates when the Distribution Pump hand switch is in the ON position.

Acid Tank Pump Run Light

The green Acid Tank Pump Run light illuminates when the Acid Tank Pump hand switch is in the ON position.

Acid Tank Low Level Light

The amber Acid Tank Low Level light illuminates when a low-level condition in the associated acid tank occurs.

PSDS Status Monitor

The PSDS Status Monitor is designed to notify personnel of the Pressurized Solution Distribution System (PSDS) operation status using a combination of visual lights and audible alarms. Monitoring logic and power are delivered from the PSDS controller through cable connections to the PSDS Status Monitor. Its compact size and optional mounting features allow it to be located near the dialysis area or nurses station. The PSDS Status Monitor is available in either an external surface mount housing or flush mount panel with a variety of monitored features that include distribution pump and distribution tank low level alarm status, acid tank levels, and alarm conditions.

1.9 Alarms

Distribution Tank High-Level Alarm

A distribution tank high-level alarm condition is created when the solution level in the distribution tank rises to engage the high-level float switch. The Mix Pump hand switch must be in the TRANSFER position and the Distribution tank high-level float switch engaged to activate the alarm. The mix pump shuts off to prevent overfilling of the distribution tank.

The high-level condition is typically observed when the solution level in the distribution tank fills during transfer from the mix tank. This is a normal occurrence and can be expected after the initial daily transfer. An audible alarm sounds to prompt the user to reposition the Mix Pump hand switch and plumbing valves once transfer is stopped.

The alarm will be discontinued when the mix pump hand switch on the PSDS is turned to the OFF position.

Distribution Tank Low-Level Alarm

A distribution tank low-level alarm condition is created when the solution level in the distribution tank drops and activates the low-level float switch. The Distribution Tank Low Level Alarm hand switch must be in the ON position and the distribution tank low-level float switch not satisfied to activate the alarm.

The low-level condition is typically observed when the solution level in the distribution tank falls during use at the end of the treatment day or prior to the initial daily batch transfer from the mix tank. An audible alarm sounds to prompt the user to recognize the low-level condition and take steps to replenish the distribution tank solution levels and/or disable the alarm. This alarm sounds when there are approximately 20 gallons of solution left in the distribution tank.

When the Distribution Tank Low Level Alarm hand switch is turned to the RESET/BYPASS position, the alarm will be discontinued.

Mix Pump Low Flow Alarm

A mix pump low flow alarm condition is created when the solution velocity in the mix pump piping drops below 2 GPM as measured at the flow switch on the output side of the mix pump. The Mix Pump hand switch must be in either the TRANSFER or MIX position and the mix pump low flow switch not satisfied to activate the alarm. The mix pump shuts off in the absence of solution to prevent pump damage.

The low flow condition is typically observed when the solution level in the mix tank empties while being transferred to the distribution tank. This is a normal occurrence and can be expected upon the initial daily transfer. An audible alarm sounds to prompt the user to reposition the Mix Pump hand switch and plumbing valves once transfer is stopped.

When the mix pump hand switch on the PSDS is turned to the OFF position the alarm will be discontinued. The Alarm Mute pushbutton switch can be pressed to silence the audible alarm prior to switch and valve repositioning.

Alarm Mute

Pressing the Alarm Mute switch will silence the audible alarm for a period of five minutes. The alarms will not sound again if the alarm condition is removed during this period with the exception of the Distribution Tank Low Level alarm that does require the resetting of the Distribution Tank Low Level hand switch.

1.10 Options

Number of Acids

Multiple product options accommodate varying customer needs for distribution loop combinations including single bicarbonate and up to three acid solutions (1 Bicarb, 0 Acid; 1 Bicarb, 1 Acid; 1 Bicarb 2 Acid; and 1 Bicarb, 3 Acid).

Wall Stations

Wall-mounted solution dispensers are capable of delivering up to three different concentrates and water as well as drain access. Process connections are configured to be compatible with dialysis manufacturer's equipment.

Bulk Acid Storage

Larger tanks (165 – 1000 gallon) are available for onsite bulk storage of acid.

Heat Disinfection (RO Water Loop):

Heat disinfection only applies to the RO water inlet piping loop that feeds the PSDS unit. The PSDS system cannot be heat disinfected.

If the PSDS unit is equipped with PVC loop flow piping, then it is **NOT** heat disinfect compatible. The PSDS unit must be equipped with PVDF loop flow piping in order to have this option. When performing RO loop piping heat disinfection, close the RO supply water shut-off valve to the PSDS unit.

WARNING:

The piping components, internal to the unit, are not heat tolerant and cannot be heat disinfected. To disinfect the PSDS system, refer to Section 3.9, PSDS Disinfection/De-Calcifying Procedure.

1.11 Specifications

1.11.1 Input Water Requirements

	<u>Minimum</u>	<u>Maximum</u>
Input flow Rate (AAMI quality water)	0.5 GPM	5.0 GPM
Temperature (water)*	10°C (50°F)	30°C (86°F)
pH (water)	3.0	11.0
Inlet Pressure	20 PSI	100 PSI
Drain Capacity	20 GPM	N/A
Drain Height	N/A	6 in. (0.15 m)

* Refer to the bicarbonate powder manufacturer's recommendation for water temperature.

1.11.2 Electrical Specifications

Voltage	115 V~
Hertz	60
Amps	20
Phase	Single

1.11.3 Environmental Requirement

	<u>Minimum</u>	<u>Maximum</u>
Ambient Temperature	4°C (39°F)	32°C (90°F)
Storage Temperature	2°C (36°F)	32°C (90°F)
Altitude	N/A	10,000 feet

NOTE:

The user/operator should recognize that moisture can be caused by condensation and is not necessarily an equipment leak. This equipment will function in the presence of condensation.

1.11.4 Dimensions/Weight

PSDS* 58"H x 110"W x 48"D

*Additional space is required for operator access and acid tank(s).

Max Operational Weight: 2050 lbs

1.12 Disposal

Disposal of this product or parts must be carried out in accordance with local disposal codes.

1.13 Electromagnetic Interference

This equipment can generate, uses and can radiate radio frequency energy and, if not installed and used in accordance with these instructions, may cause harmful interference to other devices in the vicinity. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference with other devices, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving device.
- Increase the separation between the equipment.
- Connect the equipment into an outlet on a circuit different from that to which the other device(s) is connected.
- Consult the field service technician or manufacturer for help.

1.14 Service Assistance

If service assistance is required, please take the following steps:

Consult the troubleshooting section of this manual (Chapter 5). If the problem cannot be identified and corrected by any of the procedures found in that section, then...

Contact your Facility Equipment Technician. If the technician is unable to help then...

Call the Technical Support Department at (800) 633-3080. Technicians are available for all calls between 7:00 a.m. and 7:00 p.m. CST, Monday through Friday. Technicians are also available at other times for **emergency calls only***. Product consultants will be on hand to discuss the problem with you and endeavor to rectify it over the phone. If the problem appears to be of a more serious nature, you will be given instructions regarding the action to be taken. Prior to making the phone call, you must be prepared to answer two questions:

1. What unit do you have, i.e. Pressurized Solution Distribution System
2. What is the serial number of your unit?

SERIAL NUMBER: _____

In addition, for **Non-Emergency** issues, you may e-mail technicalservice.marcor@xylem.com and a Technician will respond generally within one working day.

***Emergency:** A situation that has or could cause the shut down of a clinic during operating hours or that puts patients in increased jeopardy.

1.15 Return Material Authorization (RMA) Procedure

If you wish to return goods for, warranty evaluation and/or credit please have your original sales order, invoice and device serial number available when you call. Call us at (800) 633-3080 and request Technical Support. A representative will provide instructions and a return authorization number, which needs to be clearly written on the outside of the box used to ship your materials. All equipment must be shipped with the freight prepaid by the customer. Call our Customer Service Center with any questions or issues concerning freight claims and a representative will discuss your situation.

All materials to be returned must be rendered into a non-hazardous condition prior to shipping.

1.16 Symbols and Abbreviations

~	Volts Alternating Current
AAMI	Association for the Advancement of Medical Instrumentation
C	Celsius
CSA	Canadian Standards Association
CC	Cubic Centimeters
cfu/mL	Colony Forming Units per milliliter
EU	Endotoxin Units
F	Fahrenheit
FDA	Food and Drug Administration
GPM	Gallons Per Minute
HZ	Hertz
LBS	Pounds
LPM	Liters Per Minute
mEq	Milli equivalents
mg/L	Milligrams Per Liter
mL	Milliliters
MNPT	Male National Pipe Thread
N/A	Not Applicable
ng/mL	Nanograms Per Milliliter
NPT	National Pipe Thread
PLC	Programmable Logic Controller
P/N	Part Number
PPM	Parts Per Million
PSDS	Pressurized Solution Distribution System
PSI	Pounds Per Square Inch
PVC	Polyvinyl Chloride
QD	Quick Disconnect
RO	Reverse Osmosis
VAC or V~	Volts Alternating Current

Pressurized Solution Distribution System PSDS

CHAPTER 2: INSTALLATION GUIDELINES AND INSTRUCTIONS

This page intentionally left blank.

2.1 Installation Guidelines

This chapter provides basic guidelines for the installation of the PSDS. Each application/installation will have particular characteristics that require individual attention when planning the installation.

2.2 PSDS Installation Requirements

The Pressurized Solution Distribution System should be located close to the patient care area to minimize loop length. The equipment only footprint is 3 feet deep by nine feet wide. An additional 48" depth for operating space in front of the equipment, 24" on at least one side, and at least 24" behind are recommended. The floor area supporting the equipment must be level for proper function.

2.2.1 Electrical Requirements

115 VAC, 20 AMP, 1 Phase, 60Hz, Straight Blade Plug GFI Protected Outlet

115 VAC, 20 AMP, 1 Phase, 60Hz, GFI Protected Duplex Electrical Outlet is required for the acid distribution option.

2.2.2 Water Supply

The PSDS requires water that meets AAMI or local hemodialysis water quality standards. Additionally, it is necessary to have adequate flow rates and pressures. The minimum flow rate required for operation of the PSDS is ½ gallon per minute.

The supplied water temperature affects the bicarbonate mixing ability. Follow the bicarbonate manufacturer's recommendations for preferred water temperature.

2.2.3 Drain (Floor Sink)

The PSDS requires a drain outlet with a minimum capacity of twenty gallons per minute (20 GPM) of continuous flow and a maximum drain height of six inches.

The PSDS drain(s) should be routed separately from other equipment to a drain.

2.3 Distribution Loop Installation Requirements

Distribution loop tubing:

- Material: Pipe or Hose, ¾" I.D.
- Material: Polyethylene Tube, 5/8" O.D.

The recommended maximum number of patient stations per loop is 25 stations with a loop length not more than 600 feet. The PSDS can support a maximum of two loops.

The distribution loop starts at the distribution pump and is routed to the patient stations before returning to the distribution tank.

The distribution loop may be routed in the ceiling (or raceway) dropping down to the patient stations or under the floor, rising to the patient stations. The loop can be elevated to a maximum height of 15 feet from the pumps.

Distribution loop plumbing must be continuous (no in-line couplings in walls or ceiling). Breaks or splices in the loop should occur with tee fittings only at a patient station dispenser.

Distribution loop tubing must be installed with no kinks or pinched sections and should be installed straight (level) with appropriate restraints to prevent excessive sagging.

If the distribution loop is to pass through an unheated space, the tube/piping should be insulated to prevent freezing.

2.4 Patient Stations

Patient station installation needs will vary with the particular distribution loop material type used and the number of distribution loops. The distribution loops will utilize either ¾" ID pipe (or hosing), or 5/8" tubing, and require applicable connection hardware. Multiple distribution loops and the connections to the patient stations will require that each one is thoroughly labeled or otherwise properly identified.

The connection between the distribution loop and the patient stations should be as short as reasonably possible and must utilize FDA approved materials that are inert and compatible with disinfection methods. The organization recommends that over pressure protection (pressure regulators) is installed at each use point.

Each application/installation will have particular characteristics that require individual attention during the installation.

2.5 Installation Instructions

1. Locate the PSDS frame, mix tank, and distribution tank in a level location. The location must include a floor drain, good ventilation, a 20 amp 115 VAC 1 phase straight blade outlet, and a 20 amp 115 VAC duplex outlet if an acid distribution system is included. A water supply valve should be in reasonable proximity.
2. Clean any shipping debris from inside the mix tank and distribution tank.
3. Connect the mix tank and distribution tank assemblies to the PSDS frame.
 - a. Reference the PSDS Mix Tank To Frame Assembly and PSDS Distribution Tank To Frame Assembly illustrations in Chapter 6.
 - b. Ensure the tanks are sitting level on their stands and plumbing valves are accessible.
 - c. Connect the frame to both tanks at the matching unions and adapters. Make certain that the union O-rings are in place.
4. Assemble the drain manifold.
 - a. Reference the PSDS Drain Assembly illustration in Chapter 6.
 - b. Mount the drain manifold in a manner appropriate for the location.
 - c. The distribution loop drain connection should be constructed to allow ease of use. The bicarbonate distribution loop return will alternate routes between the drain and the distribution tank return via a three way- valve. A sample valve is included on the distribution loop return and should be positioned at an angle for gravity feed.
5. Tighten all electrical screws in the control box to 9 inch-pounds with the exception of the PLC and terminal block. The PLC terminal screws should be tightened to 5-inch pounds and the terminal block to 7 inch-pounds.

CAUTION:

Vibrations during shipping/operation may cause connections to loosen. Loose connections may not allow the PSDS to operate correctly.

6. Connect the distribution tank float switch wires to the PSDS controller.
 - a. Connect the distribution tank high-level float switch to terminal block segments COM and 4.
 - b. Connect the distribution tank low-level float switch to terminal block segments COM and 10.

7. Connect the inlet water supply.
 - a. For cold water, chemical disinfection; the unit includes two, 1 inch threaded female NPT ports for inlet and outlet RO water connections. For hot water, heat disinfection, the unit includes two $\frac{3}{4}$ inch threaded female NPT ports. Inlet and outlet connections are interchangeable at the unit. The inlet port is field connected from the loop supply feed. The outlet port is field connected to the loop return line.
8. Route the power cord to the appropriate electrical outlet.
9. Install the carbon filter for the mix and distribution tanks according to the instruction sheet in the carbon filter kit.
10. If necessary, secure the tank stands to the floor. If local codes require seismic restraints use holes in the stand's feet to secure them to the floor.

2.6 Bicarb Distribution Loop Connections

Reference the Bicarb Distribution Loop Connections illustration in Chapter 6 for visual aid.

2.6.1 Bicarb Distribution Loop Supply Connection

1. Connect the start of the bicarbonate distribution loop to the pump outlet assembly.
 - a. The installation kit includes parts to ease connection of the loop supply and not all parts of this kit may be necessary if the loop is already constructed with flexible hosing.
 - b. If a flexible connection is required use several feet of $\frac{3}{4}$ " hose, a couple of PVC nipples, hose clamps, and a $\frac{3}{4}$ " union at the necessary length to connect the distribution pumps and loop supply.

2.6.2 Bicarb Distribution Loop Return Connection

1. Connect the end of the bicarbonate distribution loop to the 3-way return valve on the distribution tank.
 - a. The installation kit includes parts to ease connection of the loop return and not all parts of this kit may be necessary if the loop is already constructed with flexible hosing.
 - b. If a flexible connection is required, use several feet of $\frac{3}{4}$ " hose, male and female quick disconnects, a PVC nipple, hose clamps, and a $\frac{3}{4}$ " union at the necessary length to connect the distribution loop return and tank assembly.

2.7 Acid Pump Controller Installation

Reference the PSDS wiring diagrams, P/N W3T577314 and W3T576135 for a visual description of wire connection locations.

NOTE: A filter and housing are provided to vent the acid tank.

2.7.1 Installation Considerations

1. Mount the Acid Pump Controller on a visible section of wall near the acid tanks and within a few feet of an 115VAC, 1 Phase, 20 Amp, GFI protected electrical outlet.
2. Connect the power cord only when all connections have been made and installation is complete.
3. Dimensions: Height – 11.5 inches; Width – 9.25 inches; Depth – 6 inches

2.7.2 PSDS Connections

1. Connect the 9 conductor cable and wiring to the PSDS Controller. Connect the fitting to the most convenient connection hole on the bottom of the controller and remove any hole plugs as necessary.
2. Connect the acid pump control wires to the PSDS controller terminal block. The number of wires is dependent on the acid controller type and not all of the following connections may be required.
 - a. Connect the white wire to terminal '24N'.
 - b. Connect the red wire to terminal '7'.
 - c. Connect the orange wire to terminal '8'.
 - d. Connect the yellow wire to terminal '9'.
3. Connect the 9-conductor cable to the Acid Pump Controller. Route the cable along the PSDS frame and present cable path.
 - a. Cut the cable to fit the length required to position the acid pump controller as close to the acid tanks as necessary.
4. Connect the acid pump control wires to the acid pump controller terminal block 2 (A or B). The number of wires is dependent on the acid controller type and not all of the following connections may be required.
 - a. Connect the white wire to terminal TB2-4.
 - b. Connect the red wire to terminal TB2-1 (Acid 1).
 - c. Connect the orange wire to terminal TB2-2 (Acid 2).
 - d. Connect the yellow wire to terminal TB2-3 (Acid 3).

2.7.3 Acid Pump Connections

1. Depending on the number of acid pumps in the particular system, connect each set of 3 acid pump power supply wires to the respective points on the Acid Pump Controller.
 - a. Connect the black wire to terminal TB1-6, TB1-7, or TB1-8 (acid 1, 2 and 3 respectively).
 - b. Connect the white wire to terminal TB1-3, TB1-4, or TB1-5 (acid 1, 2 and 3 respectively).
 - c. Connect the green wire to the ground terminal on the controller chassis.
2. Select a location for the acid pumps between the bulk tanks and the acid pump controller. Mount the acid pumps on the stands using the hardware provided. Mount the acid pump stands on the floor using the washers and hardware provided.
3. Connect the acid pump power supply wires to the acid pump motor.
 - a. Connect the two black, two white, and two green wires together.

2.7.4 Acid Tank Float Switch Connections

1. Install the acid tank low-level float switch in the normally closed position on the side of the storage tank (closed when empty).

NOTE:

The acid tank low-level float switch is normally located at the 75 to 100 gallon level. This should be discussed with the facilities Chief Technician or a representative to determine the gallon level to install the float switch.

2. Connect each pair of acid tank float switch connection wires to the PSDS Controller as shown in the system-wiring diagram in drawings P/N W3T577314 and W3T576135.
 - a. The installation kit includes cable and wire ties to connect the float switch wire leads to the PSDS Controller.

From Acid Pump Controller	To PSDS Controller
Black TB2- 5 Low Level Tank 1	Black TB1-13
Black TB2- 7 Low Level Tank 2	
Black TB2- 9 Low Level Tank 3	
Blue TB2- 6 Low Level Tank 1	Blue TB1-19
Violet TB2- 8 Low Level Tank 2	Violet TB1-20
Brown TB2- 10 Low Level Tank 3	Brown TB1-21
Red TB2-1 ADP1	TB1-7
Orange TB2-1 ADP2	TB1-8
Yellow TB2-3 ADP3	TB1-9
White TB2-4 24N	TB1-24N

2.8 Acid Distribution Loop Connections

Reference the Acid Distribution Loop Connections illustration in Chapter 6 for visual aid.

2.8.1 Acid Distribution Loop Supply Connections

1. Connect the start of the acid distribution loop to the union on the acid distribution pump outlet assembly.
 - a. The installation kit includes parts to ease connection of the loop supply and not all parts of this kit may be necessary if the loop is already constructed with flexible hosing.
 - b. If a flexible connection is required use several feet of $\frac{3}{4}$ " hose, a couple of polypropylene nipples, hose clamps, and a $\frac{3}{4}$ " union at the necessary length to connect the acid distribution pump and loop supply.

2.8.2 Acid Distribution Loop Return Connections

1. Connect the end of the acid distribution loop to the return valve on the acid distribution tank.
 - a. The installation kit includes parts to ease connection of the loop return and not all parts of this kit may be necessary if the loop is already constructed with flexible hosing.
 - b. If a flexible connection is required; use several feet of $\frac{3}{4}$ " hose, male and female quick disconnects, a polypropylene nipple, hose clamps, and a $\frac{3}{4}$ " union at the necessary length to connect the distribution loop return and tank assembly.

2.9 PSDS Status Monitor Installation

Reference the PSDS wiring diagram in Chapter 6 for a visual description of wire connection locations.

2.9.1 Installation Considerations

1. Operation of the PSDS Status Monitor is dependent on the proper connections and correct functioning of the connected devices. These devices may include the distribution tank level switches, acid tank level switches, mix pump flow switch, and the PSDS controller itself.
2. Approximate Dimensions: Height – 4”, Width – 8”, Depth – 3”
3. Mount the PSDS Status Monitor in a location with consideration for the effectiveness of sight and sound paths.
 - a. The PSDS Status Monitor is available in either a Surface Mount or Flush Mount model. The Surface Mount model mounts to the wall using slots included on the enclosure. The Flush Mount model mounts to the wall inside a 4-gang junction box.
4. The PSDS Status Monitor is directly connected to the PSDS controller through two cables. The connection cables include specific insulation color-coding requirements to correctly associate installation connections.
 - a. For cable connections use Belden brand part number 8465 (5 Conductor), and Belden brand part number 8489 (4 Conductor), or their direct equivalents.

Suggested cable vendor part numbers are supplied as follows:

<u>PSDS Cable Connection</u>	<u>Suggested Vendor Part Number</u>
5-Conductor Cable	Belden Part Number 8465
4-Conductor Cable	Belden Part Number 8489

2.9.2 Installation Procedures

1. Make wiring connections from the PSDS Status Monitor to the PSDS controller. Remove the cable jacketing as necessary.
 - a. Use the butt splice connectors on the PSDS Status Monitor to connect to the individual cables conductors. Ensure that there are no conductors visible after full tightening.

NOTE:

The 4-conductor cable will carry 24VDC signal and must be connected to the Alarm Mute switch and Distribution Tank Low Level Alarm indicator connectors.

Wiring routed from the PSDS unit to the remote monitor must meet local codes for 24V installations. Wire length should not exceed 1,000 feet.

- b. Connect the color-coded cable conductors inside the PSDS controller to their associated terminal block locations. Strip the conductor insulation for approximately ¼” to ensure a solid connection. If using the cables suggested earlier in this manual, the wire colors will match up with the wires in the monitor.

PSDS Terminal	Color	Function	Monitor Connection	Voltage
24N	White	Distribution Pump Run	J1	AC
17	Green	Distribution Pump Run	J2	AC
18	Brown	Acid Tank Low Level	J3	AC
15	Red	Distribution Tank Low Level Bypass	J4	AC
16	Black	Audio Alarm	J5	AC
Com	Black	Alarm Mute	J6	DC
14	Red	Alarm Mute	J7	DC
Com	White	Distribution Tank Low Level	J8	DC
22	Green	Distribution Tank Low Level	J9	DC

- c. Secure and cover any unused wires at the base of the removed jacketing. Not all conductors may be used depending on the monitor model and/or cable selection.

NOTES:

Pressurized Solution Distribution System PSDS

CHAPTER THREE: SYSTEM OPERATION

This page intentionally left blank.

3.1 Daily Start-Up

1. Verify the Mix Pump; Distribution Pump and Acid Pump control switches are in the OFF position; Distribution Tank Low Level in Bypass/Reset position.
2. Ensure the inlet water supply is connected to the PSDS.
3. Ensure valves V2, V3, V5, V6, V7, V10, V13, V15, V17, and V18 are closed.
4. Ensure valves V1, V4, V8 and V9 are open.
5. Ensure V14 is set to tank return.
6. Ensure the PSDS unit is plugged in to power.
7. Turn the PSDS main power switch to ON.
8. Confirm the absence of residual disinfectant prior to mixing bicarbonate solution.

WARNING:

Due to the possibility of disinfectant rebound, perform a residual disinfectant test before initiating the Bicarb Mixing Procedure.

3.2 Acid Distribution Start-Up Procedure

1. Verify sufficient acid supply for the entire treatment day.
2. Turn the appropriate Acid Tank hand switch to the ON position to start the associated acid pump.
3. Verify the respective acid tank pump run light illuminates.

WARNING:

Perform the following procedure only after the recommended disinfection and rinse procedures have been completed and a NEGATIVE residual disinfectant test is attained. Due to the possibility of disinfectant rebound, perform a residual disinfectant test before initiating the bicarbonate mixing procedure.

3.3 Mix Tank Fill Procedure

1. Ensure the availability of an inlet water supply.
2. Prepare the PSDS for mix tank fill.
 - a. Ensure valves V1, V3, V4, and V5 are open.
 - b. Ensure valves V2, V6 and V7 are closed.
3. Ensure the Main Power switch is in the ON position. Ensure the Distribution Tank Low Level Alarm switch is set to BYPASS/RESET.
4. Determine the needed amount of bicarbonate solution.
 - a. Identify how much Bicarb will be required.
 - b. Identify how much water (gallons) will be required.
5. Turn the Mix Tank Auto Fill hand switch to the MANUAL position and hold.
6. Verify water is flowing into the mix tank.
7. Set the flow meter valve, V8, to an appropriate fill rate between 0.5 and 5.0 gallons per minute. This rate will be used in combination with the Automatic Fill Timer to provide the identified water volume.
8. Release the Mix Tank Auto Fill hand switch.
9. Prepare the PSDS for mix tank fill.
 - a. Close valve V5 when tank has drain fully.
10. Verify the Auto Fill Mode Timer is correctly queued and ready for time setting.
 - a. The first identifier should be the letter "E" followed by the three numerical minute position settings and the letter "m". An illustration is provided in Chapter 6 for additional clarification.
 - b. The Auto-Fill Mode Timer may be cancelled and reset at any time during the countdown period by turning the Mix Tank Auto Fill switch to the START position.
 - c. Resetting the timer during the countdown period will result in an incorrect volume of water in the mix tank. If the timer is initiated, the fill timer will operate for the original timer setting.
11. Set the Auto Fill Mode Timer for the length of time needed in combination with the previously set fill rate to fill the mix tank to the identified water volume.

Example: 3 GPM x 20 minutes = 60 gallons
12. Turn the Mix Tank Auto Fill switch to the START position and release. The timer will blink "ON".

NOTE:

As the timer counts down the display will be reduced from a full black bar to a dashed white bar and OFF OUTPUT will be shown when finished.

13. Verify the fluid level in the mix tank once filling is complete.
 - a. If additional water is required, the Mix Tank Fill hand switch may be held in the Manual position until the needed level is reached, or the timer may be reset with a calculated fill rate and the switch set to the START position.
 - b. If the water level is greater than what is required the mix tank drain valve, V5, can be opened until the desired level is reached.
14. When tank is filled to the proper water level, close valve V3.

WARNING:

The Inlet Water valve (V3) should be closed at all times except when the mix tank is being filled with water. The Transfer Valve (V7) should remain closed at all times except when intentional transfer of solutions is to occur. If these valves are left open during treatments, and the inlet solenoid valve fails, the bicarbonate solution will be diluted, causing interruption of treatment.

3.4 Bicarb Solution Mix Procedure

Please read and follow all the bicarbonate powder manufacturers' instructions and labels regarding the preparation of bicarbonate solution before continuing with this procedure.

WARNING:

Thoroughly follow all manufacturer recommendations for mixing and testing the bicarbonate solution.

WARNING:

Valve V7 should only be OPEN during intentional transfer of fluid between mix and distribution tanks. Failure to close V7 at all other times may cause dilution or contamination of tank solutions and interruption of treatment.

NOTE:

The Mix Pump will operate for a maximum of 10 minutes. If additional mixing is required, cycle the hand switch from Mix to OFF to MIX to restart pump.

1. Ensure the mix tank has been filled to the desired solution level.
2. Ensure Valve V3 is closed.
3. Turn the Mix Pump hand switch to the MIX position and verify the water is circulating.
4. Open the mix tank lid and slowly add the appropriate amount of bicarbonate powder to the mix tank. Adding the powder to the water slowly will enable dissolving and mixing of the solution more efficiently.
5. Close the mix tank lid and allow the solution to thoroughly mix.
6. Verify the solution is completely mixed per bicarb powder manufacturer's recommendations.
 - a. A sample can be taken from the mix tank sample valve V6.
 - b. Continue mixing until the solution is completely mixed and meets all manufacturer's recommendations.

7. Turn the Mix Pump hand switch to the OFF position.
8. Bicarb Solution is ready to be transferred to the distribution tank.

3.5 Bicarb Solution Transfer Procedure

1. Ensure the bicarbonate solution has been tested and is ready for use.
2. Position the PSDS valves for distribution tank fill.
 - a. Open valve V7.
 - b. Close valve V1.
 - c. Ensure valves V4 and V9 are open.
 - d. Ensure valves V2, V3, V5, V10, V15, V17 and V18 are closed and close V13.
 - e. Ensure valve V14 is set to tank return (handle towards tank).
3. Turn the Mix Pump hand switch to the TRANSFER position and verify the solution begins transferring.
4. The solution will continue to transfer until either the distribution tank high-level float switch is engaged, or the mix pump low flow alarm is activated. Either 'a.' or 'b.' below will apply at this time.
 - a. If the Distribution Tank high level float switch is engaged (tank full), press the mute switch to silence the alarm and turn the Mix Pump Hand Switch to the OFF position and close valve V7. Any remaining bicarbonate solution remaining in the Mix Tank may be transferred later when the Distribution Tank solution is partially consumed.
 - To transfer remaining solution later, open Valve V7 and turn Mix Pump Hand Switch to the Transfer position. Complete the transfer and close valve V7.
 - b. If the Mix Tank goes empty, press the mute switch to silence the alarm and turn the Mix Pump Hand Switch to the OFF position. Close valve V7.
 - A second batch of bicarbonate solution can be mixed if required. Refer to 'Second-Batch Mix Procedure'.
5. Turn the Distribution Tank Low Level Alarm switch to the ON position.

NOTE:

The amber switch light will remain illuminated only when switch is in the Bypass/Reset position as a reminder to the operator that the Low-Level alarm is disabled.

WARNING:

Valve V7 should only be OPEN during intentional transfer of fluid between mix and distribution tanks. Failure to close V7 at other times may cause dilution or contamination of tank solutions and interruption of treatment.

3.6 Bicarb Solution Distribution Procedure

NOTE:

Ensure the bicarbonate solution has been tested and is ready for use and distribution loop ports are ready for normal operation.

1. Position the PSDS valves for distribution loop purge.
 - a. Ensure valves V11 and V13 are open.
 - b. Ensure valve V7 is closed.
 - c. Rotate valve V14 to drain (handle will point towards distribution loop).
2. Turn the Distribution Pump hand switch to the ON position and flush 10 gallons of solution into the distribution loop. This process flushes air and water to drain from loop.

CAUTION:

Monitor solution level in Distribution Tank. Do not leave unattended during this time.

3. Rotate valve V14 to the "TANK" position and continue to run pump.
4. Continue to run pump for 5 minutes to purge any remaining air and/or solution from the distribution loop and patient stations.
5. Draw a sample from Valve V17 and perform a final verification of the bicarb solution.
6. Ensure the Distribution Tank Low Level Alarm hand switch is in the ON position.
7. Set valve V13 for normal operation. Refer to either step a. or b. below:
 - a. Wall Stations WITH pressure regulators: Loop pressure adjustment is not required.
 - b. Wall Stations WITHOUT pressure regulators: Loop pressure must be set to desired loop operating pressure by adjusting valve V13. Maximum pressure must be set NOT to exceed dialysis machine manufacturer's recommendations. Throughout the day, additional adjustment may be required to maintain proper pressure.

WARNING:

Exceeding the dialysis machine manufacturer's recommended maximum pressure may cause operational or component failures. Therefore, the organization recommends the use of pressure regulation at each point of use to prevent these issues.

8. The PSDS Distribution tank and loop is now configured for normal operation for the day.

NOTE:

Monitor the Distribution tank level periodically throughout the day. An audible alarm will sound to alert the operator when the Distribution tank is low.

NOTE:

An adequate supply of bicarbonate solution must be constantly available during distribution operations. Refill the Distribution tank from the Mix tank as necessary and make additional batches of bicarb as necessary.

9. After completion of patient treatments, turn the Distribution Pump hand switch to the OFF position to discontinue distribution operations.

NOTE:

System should be drained and rinsed. PSDS should be disinfected or configured for next days' operation, depending on facility protocol.

10. All Acid Pump hand switches may be turned to OFF at this time.

3.7 Second-Batch Mix Procedure

NOTE:

After the Mix Tank has been emptied (transferred) of bicarb, an additional batch may be made.

WARNING:

Prior to transferring a second batch of bicarb, the solution must be tested and verified as correct (per bicarb manufacturer's recommendations) before beginning transfer operations.

3.7.1 Mix Tank Rinse Procedure (prior to Second Batch)

1. Open valves V2, V3, and V5.
2. Ensure valves V1 and V7 are closed.
3. Turn on the Mix Tank Auto fill switch twice, once to re-set timer and once to start the water flow. Ensure enough water is flowing through flow meter to adequately spray tank interior.
4. Allow spray to continue for 2 minutes. Stop flow by turning Mix Tank Auto Fill Switch to the Start position and release.
5. Open lid and check tank for any residual bicarb; remove any residual bicarb.
6. Allow tank to fully drain.

3.7.2 Mix Tank Fill Procedure (Second Batch)

1. Prepare the PSDS for mix tank fill.
 - a. Ensure valves V3, V4, and V5 are open; Open valve V1.
 - b. Ensure valves V6 and V7 are closed; close valve V2.
2. Determine the needed amount of bicarbonate solution.
 - a. Identify how much Bicarb will be required.
 - b. Identify how much water (gallons) will be required.
3. Turn the Mix Tank Auto Fill hand switch to the MANUAL position and hold.
4. Verify water is flowing into the mix tank.
5. Set the flow meter valve, V8, to an appropriate fill rate between 0.5 and 5.0 gallons per minute. This rate will be used in combination with the Automatic Fill Timer to provide the identified water volume.
6. Release the Mix Tank Auto Fill hand switch.
7. After tank is fully drained close valve V5.
8. Verify the Auto Fill Mode Timer is correctly queued and ready for time setting.
 - a. The first identifier should be the letter “E” followed by the three numerical minute position settings and the letter “m”. An illustration is provided in Chapter 6 for additional clarification.
9. Set the Auto Fill Mode Timer for the length of time needed in combination with the previously set fill rate to fill the mix tank to the identified water volume.

Example: 3 GPM x 20 minutes = 60 gallons
10. Turn the Mix Tank Auto Fill switch to the START position and release. The timer will blink “ON”.

NOTE:

As the timer counts down the display will be reduced from a full black bar to a dashed white bar and OFF OUTPUT will be shown when finished.

11. Verify the fluid level in the mix tank once filling is complete.
 - a. If additional water is required, the Mix Tank Fill hand switch may be held in the Manual position until the needed level is reached, or the timer may be reset with a calculated fill rate and the switch set to the START position.
 - b. If the water level is greater than what is required the mix tank drain valve, V5, can be opened until the desired level is reached.
12. Close valve V3.

WARNING:

The Inlet Water valve (V3) should be closed at all times except when the mix tank is being filled with water. The Transfer Valve (V7) should remain closed at all times except when intentional transfer of solutions is to occur. If these valves are left open during treatments, and the inlet solenoid valve fails, the bicarbonate solution will be diluted, causing interruption of treatment.

3.7.3 Bicarb Solution Mix Procedure (Second Batch)

Please read and follow all the bicarb manufacturers' instructions and labels regarding the preparation of bicarbonate solution before continuing with this procedure.

WARNING:

Thoroughly follow all manufacturer recommendations for mixing and testing the bicarbonate solution.

WARNING:

Valve V7 should only be OPEN during intentional transfer of fluid between mix and distribution tanks. Failure to close V7 at all other times may cause dilution or contamination of tank solutions and interruption of treatment.

1. Turn the Mix Pump hand switch to the MIX position and verify the water is circulating.
2. Open the mix tank lid and slowly add the appropriate amount of bicarbonate powder to the mix tank. Adding the powder to the water slowly will enable dissolving and mixing of the solution more efficiently.
3. Close the mix tank lid and allow the solution to thoroughly mix.
4. Verify the solution is completely mixed per bicarb powder manufacturer's recommendations.
 - a. A sample can be taken from the mix tank sample valve V6.
 - b. Continue mixing until the solution is completely mixed.
5. Turn the Mix Pump hand switch to the OFF position.
6. Bicarb Solution is ready to be transferred to the distribution tank as necessary.

3.7.4 Bicarbonate Transfer Procedure (Second Batch)

NOTE:

Bicarbonate solution transfer can occur when the Distribution tank level is low, and tank requires re-filling. Transfer from Mix Tank to Distribution Tank may be performed at any time during the patient treatment day.

WARNING:

Valve V7 should only be OPEN during intentional transfer of fluid between mix and distribution tanks. Failure to close V7 at other times may cause dilution or contamination of tank solutions and interruption of treatment.

1. Open valve V7; Ensure valve V9 is open; Close valve V1.
2. Turn the Mix Pump hand switch to the TRANSFER position and verify the solution begins transferring.
3. The solution will continue to transfer until either the distribution tank high-level float switch is engaged, or the mix pump low flow alarm is activated. Either 'a.' or 'b.' below will apply at this time.
 - a. If the Distribution Tank high level float switch is engaged (tank full), press the mute switch to silence the alarm and turn the Mix Pump Hand Switch to the OFF position and close valve V7. Any remaining bicarbonate solution remaining in the Mix Tank may be transferred later when the Distribution Tank solution is partially consumed.
 - To transfer remaining solution later, open Valve V7 and turn Mix Pump Hand Switch to the Transfer position. Complete the transfer and close valve V7.
 - b. If the Mix Tank goes empty, press the mute switch to silence the alarm and turn the Mix Pump Hand Switch to the OFF position and close valve V7.
 - An additional batch of bicarbonate solution can be mixed if required. Refer to 'Second-Batch Mix Procedure'.

WARNING:

Valve V7 should only be OPEN during intentional transfer of fluid between mix and distribution tanks. Failure to close V7 at other times may cause dilution or contamination of tank solutions and interruption of treatment.

3.8 End of Day - Bicarb Distribution Tank & Loop Drain and Rinse Procedure

NOTE:

After completion of patient treatments, turn the Distribution Pump Run hand switch to the OFF position to discontinue distribution operations. System should be drained and rinsed. PSDS should be disinfected or configured for next days' operation, depending on facility protocol. Acid pumps may be shut down at this time also.

WARNING:

Exceeding the dialysis machine manufacturer's recommended maximum pressure may cause operational or component failures. Therefore, the organization recommends the use of pressure regulation at each point of use to prevent these issues.

PSDS - Draining the Distribution Tank

1. Ensure the patient treatments are complete and no bicarb solution is required.
2. Turn the Acid Distribution pumps OFF.
3. Turn Distribution Pump hand switch to the OFF position.
4. Turn the Distribution Tank Low Level alarm switch to Bypass/Reset.
5. Open valve V13.
6. Rotate valve V14 to drain.
7. Open Valve V18 and drain tank.

PSDS - Draining the Mix Tank

8. Open valve V5.
9. Allow mix tank to drain completely.

PSDS - Rinsing the Mix Tank

10. Prepare the PSDS for mix tank rinse.
 - a. Open valves V1, V2, V3 and ensure valve V5 is open.
 - b. Ensure valve V7 is closed.
11. Turn the Mix Tank Auto Fill hand switch to the MANUAL position and hold.
12. Verify water is flowing into the mix tank.
13. Set the flow meter valve, V8, to the maximum available flow rate. This maximizes the spray/rinse effects inside the tank.
14. Release the Mix Tank Auto Fill hand switch and turn it to START.

NOTE:

Throughout this procedure, whenever the procedure calls for water flow to be started (i.e., "Activate the Mix Tank Auto Fill timer") it may be necessary to turn the switch to START once to re-set the timer, and to START a second time to start the water flow.

15. Allow water to flow for 1 minute.
16. Turn Mix Tank Auto Fill switch to START to stop flow; allow tank to drain fully.
17. Check for any residual bicarb solids in Mix Tank and remove if present.

18. Close valve V1 for the remaining rinses.
19. Turn Mix Tank Auto Fill switch to START to re-start flow and allow to rinse for 1 full minute again before stopping; then allow tank to drain fully.
20. Repeat the above “rinse & drain” process 3 more times to fully remove residual bicarb solution from the Mix Tank.
21. Close valves V2 and V4.
22. Turn Mix Tank Auto Fill switch to START to start flow and allow to rinse for 1 minute, then stop flow.
23. Close valve V5.
24. Open valves V7, V10 and ensure valves V18 and V9 are open.
25. Turn the Mix Tank Auto Fill hand switch to the MANUAL position and hold.
26. Verify water is flowing into the Distribution Tank.
27. Set the flow meter valve, V8, to the maximum available flow rate. This maximizes the spray/rinse effects inside the tank.
28. Release the Mix Tank Auto Fill hand switch and turn it to START.
29. Allow water to flow for 1 minute.
30. Turn Mix Tank Auto Fill switch to START to stop flow; allow tank to drain fully.
31. Check for any residual bicarb solids in Distribution Tank and remove if present.
32. Close valve V9 for the remaining rinses.
33. Turn Mix Tank Auto Fill switch to START to re-start flow and allow to rinse for 1 full minute again before stopping; then allow tank to drain fully.
34. Repeat the above 1-minute “rinse & drain” process 3 more times to fully remove residual bicarb solution from the Distribution Tank.
35. Close valve V18.
36. Open valve V9.
37. Verify valve V14 is set to drain.
38. Verify the Auto Fill Mode Timer is correctly queued and ready for time setting.
 - a. The first identifier should be the letter “E” followed by the three numerical minute position settings and the letter “m”. An illustration is provided in Chapter 6 for additional clarification.
39. Set Mix Tank Auto Fill Timer to fill Distribution Tank with 100 gallons of water.
40. Turn the Mix Tank Auto Fill switch to START, initiating water flow into Distribution Tank.

NOTE:

If the Distribution Tank will not fill beyond the level of valve V14, the water may be siphoning through loop to drain. If this is the case, it will be necessary to rotate valve V14 to the "TANK" position to fill Distribution Tank. Once the tank has been filled the V14 valve must be rotated back to the drain position.

NOTE:

During this fill process, there are no alarms or overflow protection engaged. The operator must monitor the Distribution Tank as it fills.

41. Close valve V3, V7, V9 and V10
42. Verify that valves V11 and V13 are open. If valve V14 had been turned to "TANK", rotate it back to "DRAIN" at this time.
43. Turn Distribution Pump Hand Switch to the ON position. Drain 10 gallons from the Distribution Tank.
44. Open valve V15 and drain out 1 liter to flush bicarb.
45. When Distribution Tank reaches approximately 50 gallons, sample at valve V17 for residual bicarb.
 - a. If sample is NEGATIVE, rotate valve V14 to Tank. Flush all bicarb use points in delivery loop. Rotate valve V14 to drain, and empty Distribution tank fully, then turn Distribution pump OFF.
 - b. If sample is still POSTIVE, continue to flush until Distribution Tank Low Level alarm sounds, press Mute switch and turn Distribution Tank Low Level switch to Bypass/Reset. Re-test for residual bicarb.
 - If Negative, rotate valve V14 to Tank. Flush all bicarb use points in delivery loop. Rotate valve V14 to drain, and empty Distribution tank fully, then turn Distribution pump OFF.
 - If POSITIVE, continue to rinse loop to drain. Distribution Tank may have to be re-filled and repeat the rinsing procedure to achieve negative residual test.

3.9 PSDS Disinfection/Decalcifying Procedure

WARNING:

Disinfection operations are recommended at the end of each treatment day.

WARNING:

Label all wall dispenser stations and the PSDS with appropriate warning signs, such as “Do not use, contains disinfectant”.

WARNING:

Dialysis equipment must not be connected to the PSDS during decalcifying and disinfection. Verify that all distribution points of use are disconnected.

CAUTION:

The PSDS and distribution loop must be rinsed prior to disinfection/decalcifying. Reference the PSDS Rinse procedure for instruction on system rinsing.

PSDS – Filling the Mix Tank

1. Prepare the PSDS for mix tank fill.
 - a. Open valves V1, V3, V4, and V5.
 - b. Ensure valves V2, V6 and V7 are closed.
2. Turn the Main Power switch to the ON position. Turn the Distribution Tank Low Level Alarm switch to BYPASS/RESET.
3. Turn the Mix Tank Auto Fill hand switch to the MANUAL position and hold.
4. Verify water is flowing into the mix tank.
5. Set the flow meter valve, V8, to maximum available flow rate up to 5 gallons per minute. This rate will be used in combination with the Automatic Fill Timer to provide the 60 gallons required water volume.
6. Release the Mix Tank Auto Fill hand switch.
7. Prepare the PSDS for mix tank fill.
 - a. Close valve V5 after tank is fully drained.
8. Verify the Auto Fill Mode Timer is correctly queued and ready for time setting.
 - a. The first identifier should be the letter “E” followed by the three numerical minute position settings and the letter “m”. An illustration is provided in Chapter 6 for additional clarification.

9. Set the Auto Fill Mode Timer for the length of time needed in combination with the previously set fill rate to fill the mix tank to the identified water volume.
 Example: 3 GPM x 20 minutes = 60 gallons
10. Turn the Mix Tank Auto Fill switch to the START position and release. The timer will blink "ON".

NOTE:

As the timer counts down the display will be reduced from a full black bar to a dashed white bar and OFF OUTPUT will be shown when finished.

11. Verify the fluid level in the mix tank once filling is complete.
 - a. If additional water is required, the Mix Tank Fill hand switch may be held in the Manual position until the needed level is reached, or the timer may be reset with a calculated fill rate and the switch set to the START position.
 - b. If the water level is greater than what is required the mix tank drain valve, V5, can be opened until the desired level is reached.

PSDS – Chemical Mixing and Distribution

12. Select the applicable chemical for disinfection or decalcifying. A table is listed below for reference purposes. Follow manufacturer instructions or recommendations for use and handling.

DISINFECTANT	QTY	MIX	DECALCIFIER	QTY	MIX
Bleach (Household)	9.6 cups. (2.27 LTR)	1:100	Vinegar 5%	3.0 GAL (11.35 LTR)	1:20
Peracetic Acid (Minncare HD®)	2.3 LTR	1%	Citric Acid	5 lbs	1:25
37% Formaldehyde	3.0 Gal (11.35 LTR)	2%			

WARNING:

Do not mix chemicals together.

NOTE:

The quantities listed above are for addition to 60 gallons of water to achieve the percentage or ratio of disinfectant/decalcifier in water.

NOTE:

Biosan is not recommended for use in PSDS systems due to foaming.

NOTE:

PSDS may also be disinfected with Ozone if equipped with the Ozone unit for PSDS application. Refer to PSDS Ozone Manual for Disinfection procedure.

13. Position the PSDS valves for chemical mixing operations.
 - a. Ensure valves V2, V5, and V7 are closed; close valve V3.
 - b. Ensure valves V1, and V4 are open.
14. Open the mix tank lid and gradually add the appropriate amount of chemical to the mix tank. Close the Mix Tank lid.

WARNING:

Wear gloves, eye protection and protective clothing as required by your facility. Refer to the chemical manufacturer's safety recommendations.

WARNING:

Place a label prominently on the mix tank to notify others of chemical presence: "WARNING! DECALCIFER IN USE" or "WARNING! DISINFECTANT IN USE."

15. Turn the Mix Pump hand switch to the MIX position and verify the water is circulating.
16. Allow the solution to thoroughly mix (minimum 1 minute).
17. Once the solution has been mixed, position the PSDS valves for additional spray force inside the mix tank.
 - a. Open V2.
 - b. Close valve V1.

NOTE:

Open V2 before closing V1 to avoid the risk of deadheading the pump.

18. Allow the chemical solution to mix/spray for an additional minute.
19. Turn the Mix Pump hand switch to the OFF position.
20. Position the PSDS valves for chemical solution transfer to the distribution tank.
 - a. Close valve V2.
 - b. Open valve V7.
 - c. Open valves V9 and V10.
 - d. Close valve V11.
 - e. Rotate V14 is set to tank return.
21. Turn the Mix Pump hand switch to the TRANSFER position and verify the chemical solution begins transferring.
22. Transfer 50 gallons of the chemical solution into the Distribution Tank.
23. Turn the Mix Pump hand switch to the OFF position.
24. Close valve V7.
25. Open valve V11.
26. Rotate valve V14 to Drain.
27. Turn the Distribution Pump hand switch to the ON position and drain 10 gallons.
28. Rotate valve V14 the Tank position and continue to run the Rinse Pump.

29. Draw 1 liter of solution from valve V15, V17 and V6.
30. Open valves V1 and V2; ensure valve V3 is closed.
31. Turn the Mix Pump hand switch to the Mix position.
32. Open valves V5 and V18 to drain for 5 seconds.
33. Purge any air from the distribution loop and patient stations.
34. Allow the chemical solution to circulate through the distribution loop for the appropriate contact time necessary for the particular disinfectant or decalcifier.

DISINFECTANT/DECALCIFER	CIRCULATION TIME
(1:100) Bleach	½ hour minimum
(1%) Minncare HD®	2 hours minimum
(2%) Formaldehyde	2 hours minimum
(1:20) Citric Acid	½ hour minimum

CAUTION:

Confirm the presence of the disinfectant/decalcifier at all of the bicarbonate distribution loop ports during circulation. Use appropriate test methods to verify disinfectants.

35. Turn the Mix Pump hand switch to the OFF position, Open valve V5 and allow Mix Tank to drain fully.
36. Rotate valve V14 to drain position and allow Distribution pump to completely empty the Distribution Tank.

CAUTION:

Do not leave PSDS unit unattended at this time. No alarm or pump protection is engaged at this time.

37. Turn Distribution Pump hand switch to the OFF position.

PSDS – Rinsing Chemicals from the system

38. Close valves V4 and V5
39. Ensure valves V1 and V2 are open; Open valve V3.
40. Activate the Mix Tank Auto Fill switch by turning it to the START position and allow to flow for 1 minute.
41. Turn the Mix Tank Auto Fill switch to the START position to stop flow.
42. Open valve V4 and V5 and allow to drain fully.
43. Close valves V1, V4 and V5.
44. Activate the Mix Tank Auto Fill switch by turning it to the START position and allow to flow for 1 minute.
45. Open valves V4 and V5; drain tank fully, then close valves V4 and V5.

46. Repeat the above “rinse & drain” process 3 more times to fully remove residual chemical solution from the Mix Tank.
47. Close valves V2 and open V5.
48. Start water flow for 1 minute; Stop, close valve V5.
49. Open valves V1, V2 and V4.
50. Activate the Mix Tank Auto Fill switch by turning it to the START position and fill Mix Tank to 25 gallons mark.
51. Check valve V6 for NEGATIVE residual of chemical.
 - If negative, proceed to Distribution Tank rinse procedure.
 - If POSITIVE, Drain Tank and repeat rinse steps until negative residual is achieved, then fill tank to 25 gallons.
52. Close valves V1, V2, V3 and ensure valve V5 is closed.
53. Open valve V7; Ensure valve V14 to drain. Ensure valves V4, V9 and V10 are open.
54. Turn Mix Pump hand switch to Transfer and run until pump stops and alarm sounds.
55. Press Mute switch and turn Mix Pump hand switch OFF.
56. Close valve V9.
57. Ensure valves V11 and V13 are open.
58. Turn the Distribution pump hand switch to the ON position and run until tank runs dry.

CAUTION:

Do not leave PSDS unit unattended at this time. No alarm or pump protection is engaged at this time.

59. Turn Distribution Pump hand switch to the OFF position.
60. Open valve V3.
61. Ensure valves V1, V2, V5 are closed; close valve V4.
62. Activate the Mix Tank Auto Fill switch by turning it to the START position and fill Distribution Tank to the 30-gallon mark.

NOTE:

If the Distribution Tank will not fill beyond the level of valve V14, the water may be siphoning through loop to drain. If this is the case, it will be necessary to rotate valve V14 to the “TANK” position to fill Distribution Tank. Once the tank has been filled the V14 valve must be rotated back to the drain position.

NOTE:

If valve V14 must be rotated to the TANK position, additional rinsing steps may be required.

63. Turn Mix Tank Auto Fill timer to the START position to stop water flow.

NOTE:

If valve V14 was set to TANK position in order to fill tank, valve V14 must be rotated to Drain at this time.

64. Turn Distribution Pump hand switch to the ON position.
65. Drain the Distribution Tank completely and turn the Distribution Pump OFF.
66. Open valve V9.
67. Activate the Mix Tank Auto Fill switch by turning it to the START position and fill Distribution Tank to 100 gallons.

NOTE:

If the Distribution Tank will not fill beyond the level of valve V14, the water may be siphoning through loop to drain. If this is the case, it will be necessary to rotate valve V14 to the "TANK" position to fill Distribution Tank. Once the tank has been filled the V14 valve must be rotated back to the drain position.

NOTE:

If valve V14 must be rotated to the TANK position, additional rinsing steps may be required.

CAUTION:

Do not leave PSDS unit unattended at this time. No alarm or pump protection is engaged at this time.

68. Open valve V15 and drain 2 liters to flush valve.
69. Turn the Distribution Pump hand switch to the ON position.

NOTE:

If valve V14 was set to TANK position in order to fill tank, valve V14 must be rotated to Drain at this time.

70. Open valve V15, flush out 2 liters and test for NEGATIVE residual of chemical.
71. Continue to drain tank down to 30 gallons.
72. Open valve V17, flush out 2 liters and test for NEGATIVE residual of chemical.
 - If negative, shut system down.
 - If positive, drain Distribution Tank completely and repeat filling distribution tank to 100 gallons and follow all procedural steps until negative residual test is achieved.

PSDS – System Shutdown Procedure

73. Turn Distribution Pump hand switch to the OFF position.
74. Turn Distribution Tank Low Level Alarm switch to the Bypass/Reset position.
75. Close valves V3 and V7; open valves V4, V5 and V18.
76. Turn off any acid tank pump switches.
77. Turn off the main power switch.

3.10 PSDS Status Monitor

3.10.1 Normal Operation

1. The PSDS Status Monitor is energized when the PSDS controller is connected to power and the main power switch is turned on.
2. Operation of the PSDS distribution pump is communicated with the illumination of a green Distribution Pump Run light.
3. Intentional disabling of the distribution tank low-level alarm will be observed with the illumination of the amber Distribution Tank Low Level Bypass light.
4. Any acid distribution tank low-level condition will cause the amber Acid Tank Low Level indicator to light.
5. PSDS alarms are displayed on the PSDS Status Monitor with the illumination of a red Alarm light and the sounding of an audible alarm.
6. The Alarm Mute switch allows the operator to immediately mute an audible alarm and address the situation within several minutes before re-sounding occurs. The alarm will not sound again if the alarm condition is removed during this period.

3.10.2 Alarm Condition Operation

The PSDS Status Monitor communicates alarm conditions with the illumination of the red Alarm indicator and the sounding of the audible alarm. Pressing the Alarm Mute switch on the PSDS Status Monitor will silence the audible alarm.

1. Distribution Tank Low Level Alarm

A distribution tank low-level alarm condition is created when the solution level in the distribution tank drops and activates the low-level float switch. The Distribution Tank Low Level Alarm hand switch on the PSDS must be in the ON position and the distribution tank low-level float switch not satisfied to activate the alarm.

The low-level condition is typically observed when the solution level in the distribution tank falls during use at the end of the treatment day or prior to the initial daily batch transfer from the mix tank. An audible alarm sounds to prompt the user to recognize the low-level condition and take steps to replenish the distribution tank solution levels and/or disable the alarm.

When the Distribution Tank Low Level Alarm hand switch on the PSDS is turned to the RESET/BYPASS position, the alarm will be discontinued.

2. Distribution Tank High Level Alarm

A distribution tank high-level alarm condition is created when the solution level in the distribution tank rises to engage the high-level float switch. The Mix Pump hand switch must be in the TRANSFER position and the Distribution tank high-level float switch engaged to activate the alarm. The mix pump shuts off to prevent overfilling of the distribution tank.

The high-level condition is typically observed when the solution level in the distribution tank fills during transfer from the mix tank. This is a normal occurrence and can be expected after the initial daily transfer. An audible alarm sounds to prompt the user to reposition the Mix Pump hand switch and plumbing valves once transfer is stopped.

The PSDS Status Monitor communicates a high-level condition with the illumination of the red Alarm indicator and the sounding of the audible alarm. Pressing the Alarm Mute switch on the PSDS Status Monitor will silence the audible alarm. When the mix pump hand switch on the PSDS is turned to the OFF position the alarm will be discontinued.

NOTES:

Pressurized Solution Distribution System PSDS

CHAPTER FOUR: ROUTINE MAINTENANCE

This page intentionally left blank.

4.1 Routine Maintenance of Pressurized Solution Distribution Systems

The PSDS has been designed to operate with a minimum of operator attention. Like all mechanical systems, they will operate longer with less trouble when operator maintenance is performed regularly.

See section 4.2.5 for PSDS Status Monitor Routine Maintenance.

WARNING:

Prior to performing any electrical maintenance, unplug the PSDS power cord from the wall receptacle.

NOTE:

The first part of this section provides the timelines and simple instructions for maintenance of your PSDS System. The specific instructions and procedures for the more complex maintenance items are provided in the second part of this section.

4.1.1 Daily Maintenance

Daily disinfection of the bicarbonate mixing and distribution portion of the PSDS is recommended. Refer to the PSDS Disinfection/Decalcifying procedure in chapter 3.

WARNING:

The bicarbonate mixing and distribution portion of the PSDS must be rinsed free of disinfectant, and a negative residual disinfectant test result confirmed before mixing a batch of bicarbonate solution.

4.1.2 Weekly Maintenance

Bacterial monitoring of the bicarbonate mixing and distribution portion of the PSDS is recommended on a weekly basis for at least a month after installation (or modification). Bacterial monitoring may be decreased to once a month if the weekly test results are within acceptable limits through the first month of testing.

4.1.3 Monthly Maintenance

1. Perform monthly bacterial monitoring of the bicarbonate mixing and distribution portion of the PSDS. Samples should be taken both before and after the systems have been disinfected.
2. Clean the PSDS unit frame with water to remove any acid and/or bicarbonate residue.
3. Inspect the mix and distribution tank filter for any moisture. Replace the filter if it is wet. Inspect the lines to the filter for water or any obstructions. Remove any water or obstructions found.

4.1.3.1 Leak Check

If a leak is found between threaded connectors refer to the special maintenance instructions in this chapter. As well as checking the connections for leaks, check components near or below fittings for signs of leaks.

Acid and Bicarb Distribution Loops and Patient Station Connections

1. When checking the patient stations and loop fittings for leaks, ensure that no liquid or crystal formations are present at any connection.
2. If a leak is found between the tubing and fitting cut approximately $\frac{3}{4}$ " off the end of the tubing and reconnect the fitting.

Mix Pump

1. Check the threaded connections into the large stainless steel casing of the pump and attend to any leaks to prevent damage to other components.
2. Check the seam between the motor and pump casing. If this leaks, first tighten the pump casing bolts. If the leak persists, replace the O-ring under the clamp (O-ring, P/N W2T914563).
3. Check the mechanical seal. Any leaking from the backside of the stainless pump head is likely to be a mechanical seal leak. Replace the mechanical seal in a timely manner or further damage to the motor and pump will occur (Kit, Shaft Seal, P/N W2T914211).

Bicarb Distribution and Acid Tank Pumps

1. Check the acid and bicarbonate distribution pumps for leaks. Check the fitting and tubing connections for any sign of leakage and attend to any leaks to prevent damage to other components.

4.1.4 Quarterly Maintenance

1. The bicarbonate mixing and distribution portion of the PSDS should be decalcified quarterly. Reference the PSDS Disinfection/Decalcifying procedure in Chapter 3.

NOTE:

Decalcification may be performed more frequently depending on the rate of precipitate formation and facility discretion.

2. Check the mix tank and distribution tank lid seals for degradation. Any cracks or malformations of the O-ring should be evaluated to determine if replacement is necessary.
3. Check the PSDS Status Monitor functions.
4. Verify proper operation of the PSDS check valve by following the procedure in Technote 282.

4.1.5 Semi-Annual Maintenance

WARNING:

When performing any maintenance procedure, follow all recommended safety precautions, such as de-energizing electrical circuits, personal protective gear, etc.

1. Tighten all electrical screws in the control box to 9 inch-pounds with the exception of the PLC and terminal block. The PLC terminal screws should be tightened to 5-inch pounds and terminal block to 7 inch-pounds.
2. Check all of the lights and alarms on the control panel.
 - Turn the Main Power switch to ON. Place the Mix Pump switch in the TRANSFER mode with the mix tank empty. Verify the low flow alarm sounds.
 - In the distribution mode, raise the high-level float switch in the distribution tank and verify the alarm sounds.
 - With the distribution pump running normally, verify the Distribution Pump Run indicator is illuminated.
 - Allow the distribution tank to go dry. Place the system in distribution mode. Verify the alarm sounds and the Distribution Tank Low Level indicator illuminates.
 - To verify the mute switch, push the mute button during an alarm. This will cause the audible alarm to cease, but the visual should stay on.
 - With each acid pump running, verify the related green acid pump light is illuminated.
 - Verify the acid low level indicators illuminate when the acid is below the related low level float switch.
 - Turn the Main Power switch to OFF when testing is complete.

3. Verify proper operation of the auto fill solenoid valve. See section 5.9.1 in Chapter 5, "Troubleshooting".

4.1.6 Annual Maintenance

1. Replace the mechanical seals and O-rings in the Mix Pump.
2. Inspect pumps, piping, etc. for evidence of leakage. Repair or replace fittings or components as necessary.

4.2 Specific Maintenance Procedures and Instructions

4.2.1 Bicarb System Disinfection or Decalcification Procedure

Refer to the PSDS Disinfection/Decalcifying procedure in chapter 3.

4.2.2 Recommended Bacterial Monitoring Procedure

Bacteria culture tests should be performed at least once per month, or more often as necessary to assure that the microbial count of the water used to mix bicarbonate solution complies with AAMI standards. Consult with the testing method or laboratory for required materials.

Sampling technique recommendations:

It is best to wear gloves, long sleeves and a mask to prevent contamination of the sample with skin and respiratory bacteria. Always draw the sample at the 'worst case scenario', at the end of the day just before disinfection of the system. Drawing the test sample after disinfection proves efficacy of the procedure but does not give the operator an idea if the frequency of disinfection is appropriate.

Laboratory test method:

Consult with your lab to find out what type of sterile container to use and how much water sample is required to perform the test. Make certain that your laboratory is informed of the correct procedure for performing a concentrate solution culture (e.g., a spread plate method, not a calibrated loop technique, use tryptic soy agar or equivalent, not blood agar, use one full cc, not a portion thereof). If tested incorrectly, the results may be inaccurate. The sample should be assayed within 30 minutes, or it must be refrigerated (5°C) immediately, and must be maintained in a chilled state throughout the shipping process. Shipped samples must be assayed within a 24-hour period or results may be inaccurate.

Procedure:

1. Using aseptic technique, open the sampling port and allow approximately 100 cc or more of water to flow before taking the sample.
2. After allowing the solution to flow, hold the sterile sample container under the stream.

NOTE:

Closing the sampling port before removing the container from the stream will expose the test solution to 'un-rinsed' sides of the sampling valve and may cause an erroneous culture result.

3. Fill the container with the appropriate amount of solution and remove it from the solution stream.
4. Close the sampling port.
5. Cap the specimen container immediately.
6. Label the specimen appropriately with the following:
 - Time and date sample obtained.
 - Person who obtained specimen.
 - Test to be performed - culture.
 - Sample source - bicarbonate solution, machine number or sample port area
 - Any other pertinent information or procedures required by laboratory.

4.2.3 Threaded Connection Leak Repair

NOTE:


Any leak from any part of the system should be corrected as soon as possible. Left unattended, leaks can cause secondary damage necessitating expensive repairs.

NOTE:

The use of an FDA/NSF approved silicone sealant is allowed on threaded connections. A small layer may be applied to the threads on the fitting prior to tightening.

1. Relieve any pressure in the piping/tubing system (i.e., turn off the pump, open a drain port, etc.) before performing maintenance on any threaded fitting.
2. Remove the connecting tube or disconnect the piping at the nearest union joint.
3. Remove any old tape from the threads and rewind them with 2 or 3 wraps of new PTFE tape.
4. Insert the threaded fitting into the opening and hand-tighten. Be careful to avoid cross threading. Tighten using a wrench, only, if necessary, about one-half turn.
5. Reconnect piping, pressurize, and inspect for leaks.
6. If the leak was not corrected, de-pressurize the system and tighten an additional quarter turn while being careful not to apply excessive tightening force. Pressurize the system and inspect for leaks.
7. If the leak continues, it may be necessary to remove the piping again and apply sealant or replace the threaded fitting.

4.2.4 Bicarb Mix Pump Maintenance (EBARA) (Mechanical Seal)

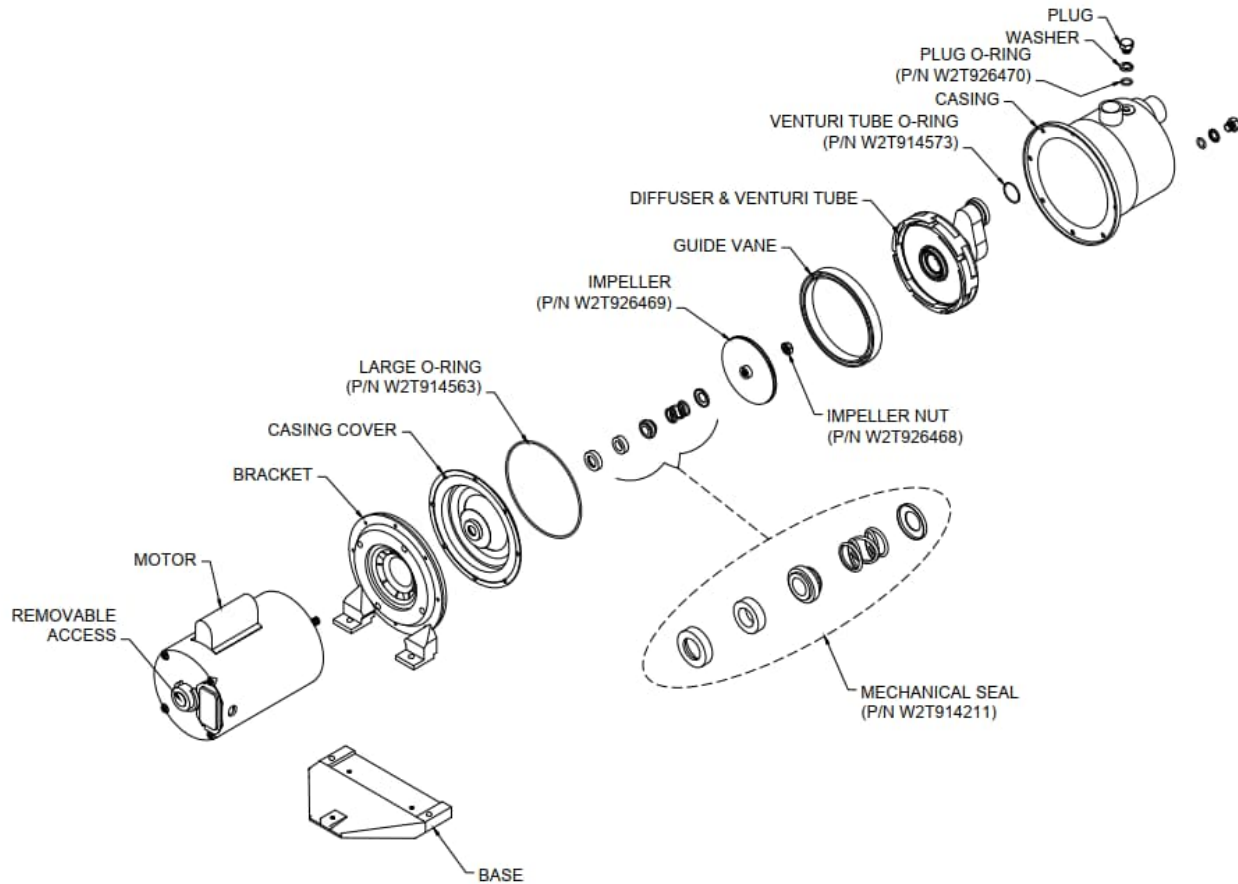

EBARA Pump Exploded View

WARNING:

Prior to performing any electrical maintenance unplug the PSDS power cord from the wall receptacle.

NOTE:

Yearly replacement of the shaft seal is recommended due the abrasive nature of bicarbonate powder.

1. Disconnect the union on the outlet and inlet piping on the pump.
2. Loosen the bolts that connect the pump casing and the motor.

WARNING:

When performing maintenance do not touch any exposed electrical connections.

3. Carefully separate the pump casing from the motor.
4. Remove the 17 mm nut holding the impeller in place. It may be necessary to grip the impeller with large channel lock pliers or a strap wrench. Use extreme caution to avoid damage to the impeller. The impeller is made of stainless steel and will tolerate 'gripping', but do not 'bite' into the metal with the pliers. If the impeller is damaged during removal, it must be replaced.
5. Remove the access cover from the rear of the motor. Use a large screwdriver to hold the shaft while removing the impeller. Unscrew the impeller from the motor shaft (counterclockwise).
6. After the impeller is removed, the shaft seal assembly will be exposed. Using two screwdrivers, pry off the rotating part of the shaft seal from the motor shaft and discard. Do not score or scratch shaft.
7. Remove and discard the stationary part of the shaft seal. It can be removed with a small screwdriver.

CAUTION:

The new seal may be damaged by the presence of body oil or lubricants on the sealing surfaces. Don't touch, mar or scratch the polished surface of either the stationary seal or the rotating seal with spring. It is advisable to handle these parts with rinsed, powder-free latex gloves. Do not allow any lubricants to come in contact with the sealing surfaces.

8. Be sure the stationary seal recess and cover plate are clean. Apply a thin coat of non-toxic silicone lube/soapy water in the surface of the central recess of the cover plate.
9. Gently push the stationary seal assembly into the recess of the cover plate. Be careful not to scratch or mar the face of the seal. Make sure the seal assembly is fully seated into the recess.
10. Clean the motor shaft to remove any debris.
11. Place the rotating assembly of the shaft seal on the motor shaft.
12. Push the assembly along the motor shaft until the 2 shaft seal faces come into contact.
13. Spin the impeller on (clockwise) until the shaft starts to turn. Place the large screwdriver in the motor access and give a quick twist clockwise to finish spinning impeller on. Ensure the impeller is securely on the shaft.
14. Secure the impeller with a strap wrench or large channel lock type pliers and install the 17 mm nut to secure the impeller in place.
15. Remove all plugs and jetted ports containing O-rings. Replace with the appropriate O-ring from the O-ring kit and reinstall into the pump housing plugs, ports and venturi ports.
16. Install the pump casing (large) O-ring onto the casing cover.
17. Replace pump casing cover. Ensure not to pinch or crush O-ring. Reinstall the eight 10 mm hex head bolts. Tighten until snug.
18. Reconnect the unions on the inlet and outlet of the pump.
19. The pump is now ready for testing to ensure there are no leaks and operation is correct.

4.2.5 PSDS Status Monitor Routine Maintenance

The PSDS Status Monitor has been designed to operate with a minimum of operator attention. Operator maintenance on the PSDS Status Monitor is limited to a daily visual and semi-annual operational check. Refer to the troubleshooting section (Chapter 5) in this manual for problems with maintenance operations.

4.2.5.1 Daily PSDS Status Monitor Maintenance

1. Observe the illumination of the Distribution Pump Run and Distribution Tank Low Level Alarm Bypass indicators during normal PSDS unit operation.
2. Observe the illumination of the Alarm indicators during the presence of any associated alarm condition.

4.2.5.2 Semi-Annual PSDS Status Monitor Maintenance

4.2.5.2.1 Distribution Tank Low Level Alarm Verification

1. Ensure the PSDS Main Power and Distribution Tank Low Level Alarm hand switches are in the ON position.
2. Drain the distribution tank to a level below the low-level float switch.
3. Verify the audible alarm sounds and the red Alarm indicator illuminates on the PSDS Status Monitor.
4. Fill the distribution tank to a level above the low-level float switch.
5. Move the Distribution Tank Low Level Alarm hand switch to the RESET/BYPASS position and then back to the ON position.
6. Verify the Distribution Tank Low Level Alarm is discontinued.
7. Turn the PSDS Main Power hand switch to the OFF position.

4.2.5.2.2 Distribution Tank High Level Alarm Verification

1. Ensure the PSDS Main Power hand switch is in the ON position and the Distribution Tank Low Level Alarm hand switch is in the RESET/BYPASS position.
2. Prepare the PSDS for solution transfer operation.
3. Turn the PSDS Mix Pump hand switch to the TRANSFER position.
4. Engage the distribution tank high-level float switch. Physically raising the float switch is an acceptable method of engaging the high-level float switch.
5. Verify the audible alarm sounds and the red Alarm indicator illuminates on the PSDS Status Monitor.
6. Turn the PSDS Mix Pump hand switch to the OFF position.
7. Verify the Distribution Tank High Level Alarm is discontinued.

4.2.5.2.3 Mix Pump Low Flow Alarm Verification

1. Ensure the PSDS Main Power hand switch is in the ON position and the Distribution Tank Low Level Alarm hand switch is in the RESET/BYPASS position.
2. Prepare the PSDS for solution mixing operation.
3. Turn the PSDS Mix Pump hand switch to the MIX position.
4. Shut the mix pump outlet valve, V1.
5. Verify the audible alarm sounds and the red Alarm indicator illuminates on the PSDS Status Monitor.
6. Turn the PSDS Mix Pump hand switch to the OFF position and open the mix pump outlet valve, V1.
7. Verify the Distribution Tank High Level Alarm is discontinued.

This page intentionally left blank.

Pressurized Solution Distribution System PSDS

CHAPTER FIVE: TROUBLESHOOTING

This page intentionally left blank.

Troubleshooting Table of Contents

SECTION	TITLE	PAGE
5.1	Alarms	5-5
5.2	Auto Fill Mode Timer	5-12
5.3	Bicarb Mixing System	5-12
5.4	Distribution Loop Problems (Acid/Bicarb)	5-15
5.5	Leaks	5-15
5.6	Main Power Switch	5-16
5.7	Programmable Logic Controller (PLC)	5-16
5.8	Inlet Solenoid Valve	5-17
5.9	PSDS Status Monitor	5-18

This page intentionally left blank.

NOTE:

The electrical schematics, flow diagrams, and other drawings included in Chapter 6, will aid in troubleshooting activities.

5.1 Alarms

5.1.1 Low Flow Alarm

1. Operation of the Low Flow Alarm is discussed in Chapter 1, General Information, and should be fully understood to ensure the troubleshooting effort is done with intent to address an unexplained activity of the alarm.
2. Ensure either the mix pump discharge valve (V1); mix tank spray head valve (V2), or distribution tank fill (V7) and distribution inlet valves (V9 and/or V10) is/are open. An open path for solution to pass through the flow switch must be available to measure flow.
3. Ensure there are no loose electrical connections between the associated controller terminal posts (TB-5 & COM), PLC terminal posts (I10 & 24VDC) and mix pump flow switch. Repair any loose electrical connections as necessary.
4. Verify the Low Flow alarm condition is registered on the PLC.
 - a. The alarm condition will be registered in the PLC with the non-illumination of the associated PLC Input LED #10. If the alarm condition is not being registered in the PLC, there may be loose wire connections or failed components including the PLC. Repair any loose electrical connections as necessary and continue to investigate possible component failure if the problem remains unresolved.
5. Verify the Low Flow switch is functional.
 - a. Ensure the Main Power hand switch is in the OFF position.
 - b. Disconnect the flow switch from its wiring and verify an infinite resistance across the flow switch leads. If an infinite resistance measurement is not observed the PLC may have failed and require replacement.
 - c. If a near zero resistance measurement is observed the switch may have been stuck or failed and require replacement.

5.1.2 Distribution Tank High-Level Alarm

1. Operation of the Distribution Tank High-Level Alarm is discussed in Chapter 1, General Information, and should be fully understood to ensure the troubleshooting effort is done with intent to address an unexplained activity of the alarm.
2. Ensure the distribution tank high-level float switch is installed in the normally closed (NC) position and has not been stuck or broken.
3. Ensure there are no loose electrical connections between the associated controller terminal posts (TB-4 & COM), PLC terminal posts (I9 & 24VDC COM) and distribution tank high-level float switch. Repair any loose electrical connections as necessary.
4. Verify the Distribution Tank High-Level alarm condition is registered on the PLC.
 - a. The alarm condition will be registered in the PLC with the non-illumination of the associated PLC Input LED #9. If the alarm condition is not being registered in the PLC, there may be loose wire connections or failed components including the PLC. Repair any loose electrical connections as necessary and continue to investigate possible component failure if the problem remains unresolved.
5. Verify the distribution tank high-level float switch is functional.
 - a. Disconnect the float switch from its wiring and verify continuity across the float switch leads with the switch in the normally closed (NC) position. If continuity is not observed, the switch may have failed and require replacement.
 - b. Verify an infinite resistance across the float switch leads with the switch in the normally open (NO) position. If a near zero resistance measure is observed, the switch may have failed and require replacement.
 - c. If correct operation of the float switch is observed, the PLC may have failed and require replacement.

5.1.3 Distribution Tank Low-Level Alarm

1. Operation of the Distribution Tank Low-Level Alarm is discussed in Chapter 1, General Information, and should be fully understood to ensure the troubleshooting effort is done with intent to address an unexplained activity of the alarm.
2. Ensure the distribution tank low-level float switch is installed in the normally open (NO) position and has not been stuck or broken.
3. Ensure there are no loose electrical connections between the associated controller terminal posts (TB-10 & COM), PLC terminal posts (I13 & 24VDC COM), and distribution tank low-level float switch. Repair any loose electrical connections as necessary.

4. Verify the Distribution Tank Low-Level alarm condition is registered on the PLC.
 - a. The alarm condition will be registered in the PLC with the non-illumination of the associated PLC Input LED #13. If the alarm condition is not being registered in the PLC, there may be loose wire connections or failed components including the PLC. Repair any loose electrical connections as necessary and continue to investigate possible component failure if the problem remains unresolved.
5. Verify the distribution tank low-level float switch is functional.
 - a. Disconnect the float switch from its wiring and verify an infinite resistance across the float switch leads with the switch in the normally open (NO) position. If a near zero resistance measure is observed, the switch may have failed and require replacement.
 - b. Verify continuity across the float switch leads with the switch in the normally closed (NC) position. If continuity is not observed, the switch may have failed and require replacement.
 - c. If correct operation of the float switch is observed, the PLC may have failed and require replacement.
6. Verify the Distribution Tank Low Level Alarm hand switch is operational.
 - a. Verify switch function with a continuity check across the switch terminals in both the normally open and normally closed switch positions. Continuity should only be observed in the normally closed position. Replace the hand switch if necessary.

5.1.4 Audible Alarm

1. Operation of the Audible Alarm is discussed in Chapter 1, General Information, and should be fully understood to ensure the troubleshooting effort is done with intent to address an unexplained activity of the alarm.
2. Ensure the Alarm Mute is not temporarily silencing the audible alarm.
 - a. Move the Main Power hand switch to the OFF position and then to the ON position. If the audible alarm does not sound upon turning the switch to the ON position, there may be loose wire connections or failed components including the Sonalert. Repair any loose electrical connections as necessary and continue to investigate possible component failure if the problem remains unresolved.

3. Verify the known alarm condition is registered on the PLC.
 - a. The alarm condition will be registered in the PLC with the non-illumination of the associated PLC Input LED# (9, 10, or 13). If the alarm condition is not being registered in the PLC, there may be loose wire connections or failed components including the PLC.
 - b. Ensure there are no loose electrical connections between the associated controller terminal posts, PLC terminal posts, and Sonalert. Repair any loose electrical connections as necessary and continue to investigate possible component failure if the problem remains unresolved.
4. Verify power is being supplied to the Sonalert.
 - a. The Sonalert's PLC Output LED # 8 should be lit in the presence of the known alarm condition. If the Output LED #8 is off, the PLC may be defective and need replacement.
5. Verify the Sonalert is functioning.
 - a. Verify continuity between PSDS controller terminal post 24N and the negative terminal (-) on the Sonalert. Repair and/or reconnect any loose wires as necessary.
 - b. Verify continuity between PLC Output LED #8, PSDS controller terminal post 16, and the positive terminal (+) on the Sonalert. Repair and/or reconnect any loose wires as necessary.
 - c. If the Sonalert's PLC Output LED #8 is lit and continuity has been verified, the Sonalert may be defective and need replacement.

5.1.5 Distribution Pump Run Light

1. Ensure the distribution pump fuse (5 Amp) is not blown and there are no loose electrical connections between the associated controller terminal posts, control relay (CR1), hand switch (HS-3), and light (LT1). Replace the fuse and/or repair any loose electrical connections as necessary.
2. Verify the Distribution Pump Run light is receiving power.
 - a. Place the Main Power and Distribution Pump hand switches in the ON position.
 - b. Verify an 115VAC measurement (approximate) across the light (LT1) terminals. If the voltage is verified and the light is not lit, it may have failed and require replacement. If the voltage is not verified there may be loose wire connections or failed components. Repair any loose electrical connections as necessary and continue to investigate possible component failure if the problem remains unresolved.
3. Verify the distribution pump control relay (CR1) is functional.
 - a. Place the Main Power switch in the ON position.

- b. Toggle the Distribution Pump hand switch several times while observing the control relay (CR1) for movement of the relay contact arm. With the distribution pump hand switch in the ON position, verify voltage across the control relay (CR1) terminals 13 and 14. If the relay arm is not moving with switch activation and/or voltage is confirmed across terminals 13 and 14, the relay may have failed and require replacement.
4. Verify the Distribution Pump hand switch is operational.
 - a. Verify switch function with a continuity check across the switch terminals in both the normally open and normally closed switch positions. Continuity should only be observed in the normally closed position. Replace the hand switch if necessary.

5.1.6 Distribution Tank Low Level Light

1. Ensure there are no loose electrical connections between the associated controller terminal posts, hand switch (HS-5), PLC, and light (LT2). Repair any loose electrical connections as necessary.
2. Verify the Distribution Tank Low Level light is receiving power.
 - a. Place the Main Power and Distribution Tank Low Level hand switches in the ON position.
 - b. Create a distribution tank low-level alarm condition or disconnect the low-level float switch wires at the tank and connect the two ends to provide a circuit jumper.
 - c. Verify a 24VAC measurement (approximate) across the light (LT2) terminals. If the voltage is verified and the light is not lit, it may have failed and require replacement. If the voltage is not verified there may be loose wire connections or failed components. Repair any loose electrical connections as necessary and continue to investigate possible component failure if the problem remains unresolved.
 - d. If applicable, reconnect the Distribution Tank low-level float switch wires once verification operations have been completed.
3. Verify the Distribution Tank Low Level hand switch is operational.
 - a. Verify switch function with a continuity check across the switch terminals in both the normally open and normally closed switch positions. Continuity should only be observed in the normally closed position. Replace the hand switch if necessary.
4. Verify the distribution tank low-level alarm is functional. Reference the Distribution Tank Low Level Alarm section of this chapter for additional troubleshooting instruction.

5.1.7 Acid Tank Pump Run Light(s)

1. Ensure the acid pump fuse (2.5 Amp) in the Acid Pump Controller is not blown and there are no loose electrical connections between the associated controller terminal posts, control relay (CR5, CR6, or CR7), hand switch (HS-6, HS-7, or HS-8), and light (LT3, LT5, or LT7). Replace the fuse and/or repair any loose electrical connections as necessary.
2. Verify the Acid Tank Pump Run light is receiving power.
 - a. Place the Main Power and Acid Tank hand switch (HS-6, HS-7, or HS-8) in the ON position.
 - b. Verify a 24VAC measurement (approximate) across the light (LT3, LT5, or LT7) terminals. If the voltage is verified and the light is not lit it may have failed and require replacement. If the voltage is not verified, there may be loose wire connections or failed components. Repair any loose electrical connections as necessary and continue to investigate possible component failure if the problem remains unresolved.
3. Verify the acid pump control relay (CR5, CR6, or CR7) is functional.
 - a. Place the Main Power switch and the Acid Tank hand switch (HS-6, HS-7, or HS-8) in the ON position while observing the associated control relay for movement of the relay contact arm. Verify voltage across the control relay (CR5, CR6, or CR7) terminals 13 and 14. If the relay arm is not moving with switch activation and/or voltage is confirmed across terminals 13 and 14, the relay may have failed and require replacement.
4. Verify the Acid Tank hand switch is operational.
 - a. Verify switch function with a continuity check across the switch terminals in both the normally open and normally closed switch positions. Continuity should only be observed in the normally closed position. Replace the hand switch if necessary.

5.2 Auto-Fill Mode Timer

5.2.1 Auto Fill Mode Timer Display Window (LCD) Is Blank

1. Operation of the Auto-Fill Mode Timer is discussed in Chapter 1, General Information, and should be fully understood to ensure the troubleshooting effort is done with intent to address an unexplained activity of the alarm.
2. Ensure there are no loose electrical connections between the associated controller terminal posts (6, 13, 24N, & COM), PLC terminal posts (I8, Q6, Q7, Input & Output Common), and the timer relay base (Posts 2, 3, 6, 7, 9, 10, & 11). Repair any loose electrical connections as necessary.
3. Ensure the power cord is connected to an 115VAC power and the main power switch is turned to the ON position.
4. Verify the Auto Fill Timer is receiving power.
 - a. Place the Main Power hand switch in the ON position.
 - b. Verify a 24VAC measurement (approximate) across terminals #2 and #10 of the timer relay base. If the voltage is verified, the Auto Fill Timer may be defective and require replacement. If the voltage is not verified, there may be loose wire connections or failed components. Repair any loose electrical connections as necessary and continue to investigate possible component failure if the problem remains unresolved.

5.2.2 No Water Flow (With Timer Operating Correctly)

1. Ensure the water supply valve, flow meter needle valve (V8), and inlet water valves (V3 and V1 or V2) are open.
2. If the valve positioning is correctly verified, the Inlet Solenoid may be defective and require replacement. Reference the Inlet Solenoid troubleshooting section in this chapter for additional verification operations.

5.3 Bicarb Mixing System

5.3.1 Mix Pump Motor Does Not Start

1. Ensure the circuit breaker, located on the left side panel of the controller, is in the ON position.
2. Ensure there are no loose electrical connections between the associated controller terminal posts (H, 1, 2, N, & 24N), PLC terminal posts (Q12 & 24VAC COM), motor contactor relay (MCR1), circuit breaker (CB-1), and motor contactor (C-1). Repair any loose electrical connections as necessary.
3. Ensure the power cord is connected to 115VAC power and the main power switch is in the ON position.
4. If the Mix Pump hand switch is in the Mix position and the motor will not start, ensure the low flow alarm is not being activated.
5. Verify the Low Flow alarm is functional. Reference the Low Flow Alarm section of this chapter for additional troubleshooting instruction.
6. If the Mix Pump hand switch is in the Transfer position and the motor will not start,

- ensure that both the low flow alarm and distribution tank high-level switch alarm are not being activated.
- a. Verify the Low Flow and Distribution Tank High-Level alarms are functional. Reference the Low Flow Alarm and Distribution Tank High-Level Alarm sections of this chapter for additional troubleshooting instruction.
7. Verify the mix pump signal is registered on the PLC.
 - a. The mix pump signal will be registered in the PLC with the illumination of PLC Input LED #2 (Mix Pump 'ON' Signal) or PLC Output LED #12 (Mix Pump 'RUN' Signal). If the mix pump signal is not being registered in the PLC, there may be loose wire connections or failed components including the PLC. Repair any loose electrical connections as necessary and continue to investigate possible component failure if the problem remains unresolved.
 8. Verify the Mix Pump hand switch is operational.
 - a. Verify switch function with a continuity check across the switch terminals in both the normally open and normally closed switch positions. Continuity should only be observed in the normally closed position. Replace the hand switch if necessary.
 9. Verify the motor contactor relay (MCR1) is functional.
 - a. Verify an 115VAC voltage (approximate) measurement across the motor contactor relay terminals '4' and '7'. If the voltage is not verified, the motor contactor relay may have failed and require replacement.
 - b. Verify a 24VAC voltage (approximate) measurement across the motor contactor relay terminals '6' and '9'. If the voltage is not verified, the motor contactor relay may have failed and require replacement.
 10. Verify the motor contactor (C-1) is functional.

WARNING:

The motor contactor contains exposed metal parts, which normally carry high voltage electrical power (115 VAC). Do not touch the contactor except as required using appropriate insulated test meter leads.

- a. Toggle the Mix Pump switch a few times while looking for movement of the C-1 relay arm. If no movement is observed, the motor contactor may have failed and require replacement.
- b. Verify a 24VAC voltage (approximate) measurement across the PLC terminal 'Q12' and controller terminal '24N'. If the voltage is not confirmed, the contactor may have failed and require replacement.
- c. Verify an 115VAC voltage (approximate) measurement across controller terminal '2' and motor contactor terminal '5' while the Mix Pump hand switch is in either the Mix or Transfer positions. If no voltage is confirmed at C1-5, make sure the circuit breaker is in the ON position. If the circuit breaker is in the ON position and no voltage is confirmed at controller terminal '2', the motor contactor may have failed and require replacement.

11. Verify an 115VAC voltage (approximate) measurement across the controller terminal blocks '2' and 'N' and at the comparable motor connections. Repair any loose connections as necessary.
 - a. If the voltage is confirmed at the motor and it does not run, it may have failed and require replacement.

5.3.2 Mix Pump Leaks

1. If the mix pump is leaking at the inlet or outlet piping, disconnect the plumbing at the unions and replace the PTFE tape on the pipe threads. Reconnect fittings and check for leaks during operation. Reference the Leak Repair sections in Chapter 4 for additional repair instruction.
2. If the pump housing (stainless steel cover) is leaking, attempt to stop the leak by tightening the pump casing bolts. Additional replacement of the pump casing O-ring may be required to stop this leak.

NOTE:

A leak in the pump housing seal is usually seen as a small puddle of fluid directly below the motor. Bicarbonate solution is corrosive to some metals, and if left unattended, will cause damage to the pump backing plate (or 'motor stool'). Quick action to repair the leak can help reduce maintenance and repair expense.

3. If the pump shaft seal (mechanical seal) is leaking, the shaft seal must be replaced. Reference the Mix Pump Maintenance section in Chapter 4 for instruction.

CAUTION:

A leak in this area can cause damage to the motor, and/or motor stool, so it is best to repair it immediately.

5.4 Distribution Loop Problems (Acid/Bicarb)

5.4.1 Dialysate Conductivity

WARNING:

Distribution loop problems may result in improper quantities of concentrate fluids delivered to the dialysis machines, which can result in low or high conductivity dialysate. If the PSDS is found to be functioning properly, refer to the dialysis equipment for conductivity issues.

1. Ensure the correct acid concentrate is in use.
2. Ensure the bicarbonate concentrate mixture is correct.
 - a. Reference the Bicarb Solution Mix Procedure in Chapter 3 for mixing operations.
3. Ensure the flow and pressure is sufficient for use.
 - a. Verify the distribution pump output valve is opened enough to allow sufficient flow.
 - b. Verify there are no kinks or air leaks in the tubing or fittings. Inspect and repair as necessary.

5.4.2 Fluid Leaks

1. Ensure the fittings and/or O-rings are not damaged. Repair or replace as necessary.
2. Replace scratched or scored tubing as necessary.

5.5 Leaks

1. A leak from any part of the system should be corrected as soon as possible. Left unattended, leaks can cause secondary damage necessitating expensive repairs. Reference the Leak Repair sections in Chapter 4 for additional instruction.

5.6 Main Power Switch

5.6.1 No Power

1. Ensure the power cord is connected to an 115VAC power supply.
2. Ensure the circuit breaker, located on the left side panel of the controller, is in the ON position.
3. Ensure the primary fuse (1 Amp) and secondary fuse (1.5 Amp) are not blown and there are no loose electrical connections between the associated controller terminal posts, circuit breaker (CB-1), motor contactor (C-1), motor contactor relay (MCR1), and hand switch (HS-1). Replace the fuse and/or repair any loose electrical connections as necessary.
4. Verify the Main Power hand switch is operational.
 - a. Verify switch function with a continuity check across the switch terminals in both the normally open and normally closed switch positions. Continuity should only be observed in the normally closed position. Replace the hand switch if necessary.
5. Verify the transformer is functional.
 - b. Verify the 115VAC primary voltage (approximate) measurement across the transformer end of the primary fuse (1 Amp) and controller terminal 'N'.
 - c. Verify the 24VAC secondary voltage (approximate) measurement across the transformer end of the secondary fuse (1.5 Amp) and controller terminal '24N'.
 - d. If the primary voltage is verified and the secondary voltage is not, the transformer may have failed and require replacement.

5.7 Programmable Logic Controller (PLC)

1. Ensure the main power switch is in the ON position and all three PLC 'PWR', 'OK' and 'RUN' LED's are illuminated. If these LED's are not illuminated the PLC may have failed and require replacement.
 - a. Ensure the Run/Program switch on the PLC is in the Run mode. If the PLC 'PWR' and 'OK' LED's are both lit with the 'RUN' LED off, the PLC may have failed and require replacement.
 - b. If the PLC 'PWR' LED is lit with the 'OK' and 'RUN' LED's off, the PLC may have failed and require replacement.

NOTE:

Other LED's may also be lit depending upon the status of the hand switches and float switches.

2. If the voltage is not verified, the PLC may have failed and require replacement.
3. Verify continuity between PLC terminal 'N' and any controller terminal 'N'. If the continuity cannot be verified, there may be loose wire connections. Repair any loose wire connections as necessary.

4. Verify the transformer is functional.
 - a. Verify the 115VAC primary voltage (approximate) measurement across the transformer end of the primary fuse (1 Amp) and controller terminal 'N'.
 - b. Verify the 24VAC secondary voltage (approximate) measurement across the transformer end of the secondary fuse (1.5 Amp) and controller terminal '24N'.
 - c. If the primary voltage is verified and the secondary voltage is not, there may be loose wire connections, or the transformer may have failed and require replacement. Repair any loose wire connections as necessary.
5. Verify the motor contactor relay (MCR1) is functional.
 - a. Toggle the Main Power hand switch a few times while looking for movement of the contact arm. If no movement is observed, the motor contactor relay may have failed and require replacement.

5.8 Inlet Solenoid Valve

5.8.1 Solenoid Valve Fails To Open When Auto Fill Timer Is Activated

1. Move the Auto Fill hand switch to the Manual position and hold it for 5 seconds. Verify the solenoid valve opens and water begins flowing through the solenoid. If the solenoid does not open several tests should be done to identify the possible operation failure mode.
2. Verify the Auto Fill hand switch is operational.
 - a. Verify switch function with a continuity check across the switch terminals in both the normally open and normally closed switch positions. Continuity should only be observed in the normally closed position. Replace the hand switch if necessary.
3. Verify a 24VAC voltage measurement across the controller terminals '13' and '24N'. If the voltage is not verified, refer to the Main Power section of this chapter for additional troubleshooting instruction.
4. Hold the Auto Fill switch to the Manual position and verify a 24VAC voltage measurement across the controller terminals '6' and '24N'. If the voltage is not verified, the solenoid may have failed and require replacement.
5. Check continuity between Auto Fill Timer base terminal #9 and controller terminal '6'. If the continuity cannot be verified, there may be a loose wire connection. Repair any loose wires as necessary.

6. Hold the Auto Fill hand switch to the MANUAL position and verify continuity between the PLC 24VDC COM and PLC Input LED #8. If the continuity cannot be verified, the PLC may have failed and require replacement.
7. Turn the Auto Fill hand switch to the MANUAL position several times and verify PLC Input LED #8 lights each time the switch is turned to the Start position. If the LED illumination cannot be verified, the electrical connections may be loose and need to be repaired.
8. Turn the Auto Fill hand switch to the START position several times and verify PLC Output LED #6 or #7 light momentarily each time Input LED #8 lights. If the LED illumination cannot be verified, the PLC may have failed and require replacement.
9. Turn the Auto Fill hand switch to the START position several times and verify the Auto-Fill Timer window displays the flashing word 'ON' when PLC Output LED #6 lights. Verify the Auto-Fill Timer window displays the word 'OFF' when PLC Output LED #7 lights. If the combination of the lights and window display cannot be verified, there may be loose wire connections, the Auto Fill Timer settings may be incorrect, or the Auto Fill Timer may have failed and require replacement. Repair any loose wire connections and/or reset the Auto Fill Timer as required. Reference the Auto Fill Timer procedure.

5.8.2 Solenoid Valve Fails To Close

1. If the solenoid valve fails to close upon completion of the filling operations, it may be stuck open. Turn the Auto Fill hand switch to the Manual position several times to clear an obstruction or free a sticky valve. Replace the solenoid valve if its failure to close upon completion of the filling operations is repeated.

5.9 PSDS Status Monitor

Reference the PSDS wiring diagram in Chapter 6 for a visual description of wire connection locations.

5.9.1 Failed Indicator Lights

1. Connections
2. Verify the indicator wire leads inside the PSDS Status Monitor are correctly and securely attached to the cable connection.
3. Verify the wire lead connections inside the PSDS controller are correctly and securely attached.
4. Restore connections as required.
5. Indicators
 - a. Verify the indicator lamps have not burnt out with a continuity check across the terminal connections.
 - b. Replace indicator as necessary.

5.9.2 Failed Audible Alarm

1. Connections
 - a. Verify the Sonalert's wire leads inside the PSDS Status Monitor are correctly and securely attached to the cable connection.
 - b. Verify the Sonalert's wire lead connections inside the PSDS controller are correctly and securely attached.
 - c. Restore connections as required.
2. Sonalert
 - a. Verify the Sonalert is operative with the application of a 24 VAC power source across the terminal connections.
 - b. Replace Sonalert as necessary.

5.9.3 Failed Mute Switch

1. Connections
 - a. Verify the mute switch wire leads inside the PSDS Status Monitor are correctly and securely attached to the cable connection.
 - b. Verify the mute switch wire lead connections inside the PSDS controller are correctly and securely attached.
 - c. Restore connections as required.
2. Mute Switch
 - a. Verify the 'momentary on' mute switch is operative with a continuity check while pressing the mute switch.
 - b. Replace switch as necessary.

Pressurized Solution Distribution System PSDS

CHAPTER SIX: DRAWINGS

This page intentionally left blank.

Illustration Notes:

1. This manual includes drawings that cover available components and configurations. Some drawings may not be applicable to your unit.
2. The PSDS controller electrical schematic and system wiring diagram are included at the end of this manual in a larger fold out format to aid readability. The PSDS controller electrical schematic includes notes specifying the differences of the multiple acid configurations. The PSDS wiring diagram includes connections for the multiple acid configurations.

List of Illustrations:

6-1	PSDS Assembly
6-2	PSDS Valve Numbers and Locations
6-3	PSDS Bicarb Loop Flow Schematic
6-4	PSDS Acid Loop Flow Schematic
6-5	PSDS Controller Display Panel – Front View
6-6	PSDS Controller Display Panel – Rear View
6-7	PSDS Controller Components – Inside
6-8	GE VersaMax PLC
6-9	Electrical Component – Control Relay Base and Motor Contact
6-10	Electrical Component – Auto Fill Timer
6-11	Mix Pump/Motor Detail
6-12	Mix Pump Electrical Configuration
6-13	PSDS Mix Tank
6-14	PSDS Distribution Tank
6-15	PSDS Mix Tank to Frame Assembly
6-16	PSDS Distribution Tank to Frame Assembly
6-17	PSDS Drain Assembly
6-18	PSDS Bicarb Loop Connections
6-19	PSDS Acid Distribution Loop Connections
6-20	PSDS Acid Tank Connections
6-21	PSDS Status Monitor, Surface Mount Model
6-22	PSDS Status Monitor, Flush Mount Model
6-23	PSDS Wiring Diagrams (Located at the end of manual)

Illustration 6-1: PSDS Assembly

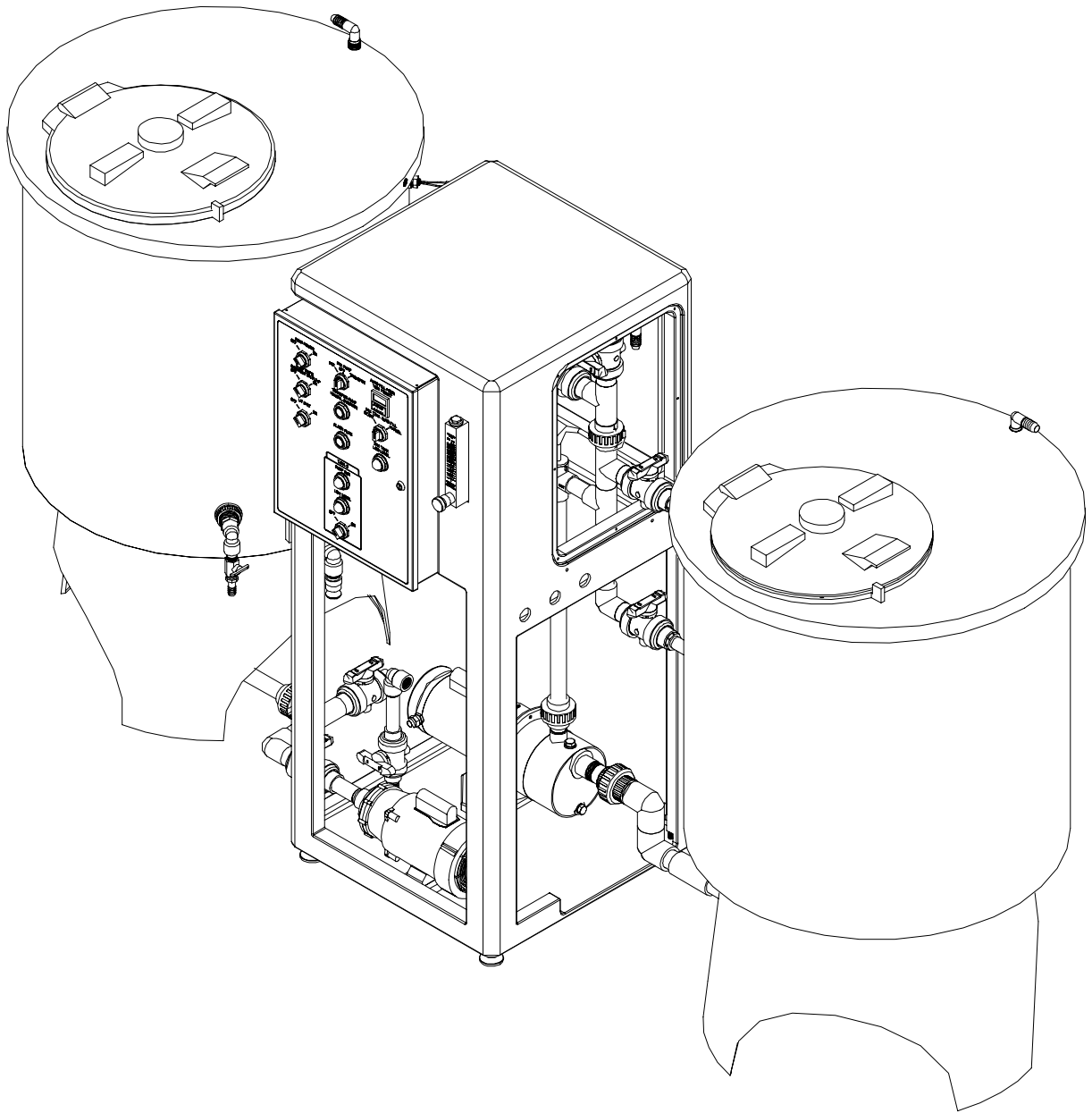


Illustration 6-2: PSDS Assembly Valve Labeling



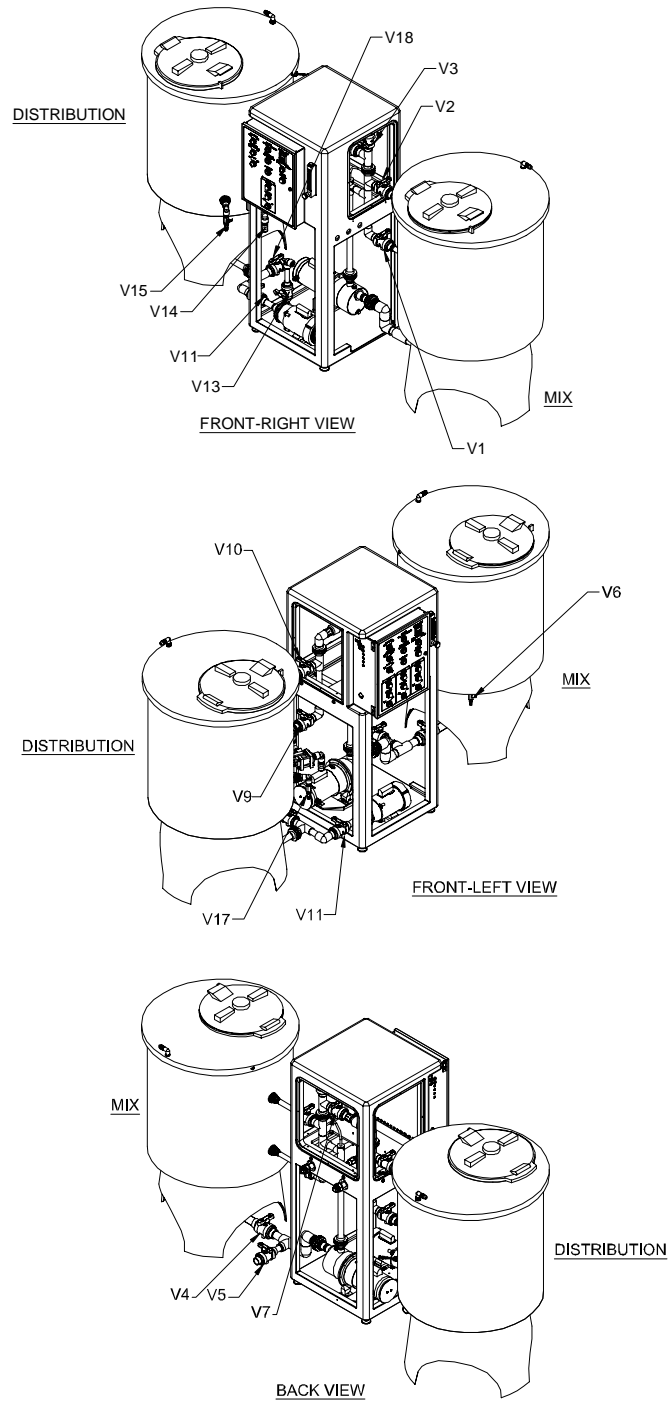


Illustration 6 - 3: PSDS Bicarb Loop Flow Schematic

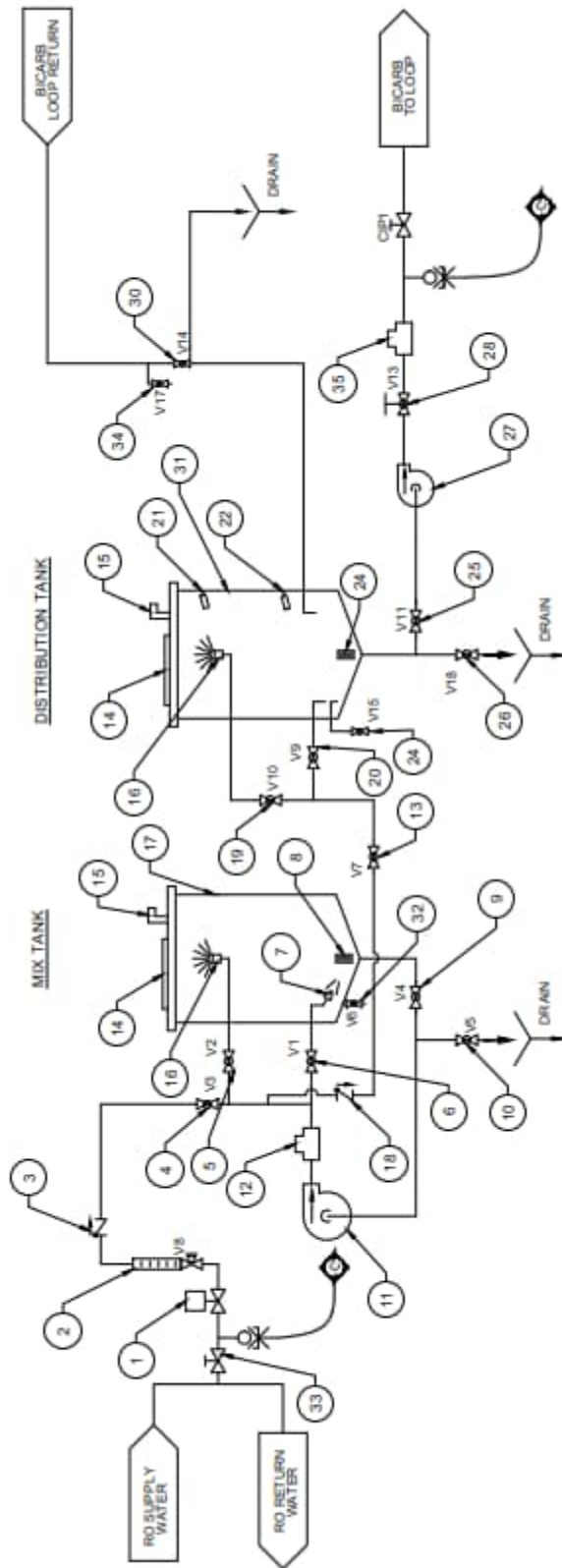


ITEM	DESCRIPTION
17	MIX TANK
16	SPRAY HEAD
15	VENT FILTER
14	TANK LID AND SEAL
13	DISTRIBUTION TANK FILL VALVE (V7)
12	FLOW SWITCH
11	MIX PUMP
10	MIX TANK DRAIN VALVE (V5)
9	MIX PUMP SUCTION VALVE (V4)
8	MIX TANK DRAIN
7	MIX EXTRACTOR
6	MIX PUMP DISCHARGE VALVE (V1)
5	MIX TANK SPRAY HEAD SUPPLY VALVE (V2)
4	INLET WATER VALVE (V3)
3	CHECK VALVE
2	ADJUSTABLE INLET FLOW METER (V6)
1	INLET SOLENOID
35	DISTRIBUTION PUMP FLOW SENSOR
34	DISTRIBUTION LOOP SAMPLE VALVE (V17)
33	SUPPLY WATER SHUTOFF VALVE
32	MIX TANK SAMPLE VALVE (V6)
31	DISTRIBUTION TANK
30	DISTRIBUTION RETURN VALVE (V14) 3 WAY
28	DISTRIBUTION FLOW CONTROL VALVE (V13)
27	DISTRIBUTION PUMP
26	DISTRIBUTION TANK DRAIN VALVE (V18)
25	DISTRIBUTION PUMP SUCTION VALVE (V11)
24	BICARBONATE SUPPLY VALVE (V15)
23	DISTRIBUTION TANK DRAIN
22	DISTRIBUTION TANK LOW LEVEL SWITCH (NO)
21	DISTRIBUTION TANK HIGH LEVEL SWITCH (NO)
20	DISTRIBUTION INLET VALVE (V9)
19	DISTRIBUTION TANK SPRAY HEAD SUPPLY VALVE (V10)
18	DISTRIBUTION FILL CHECK VALVE

Illustration 6-4: PSDS Acid Loop Flow Schematic

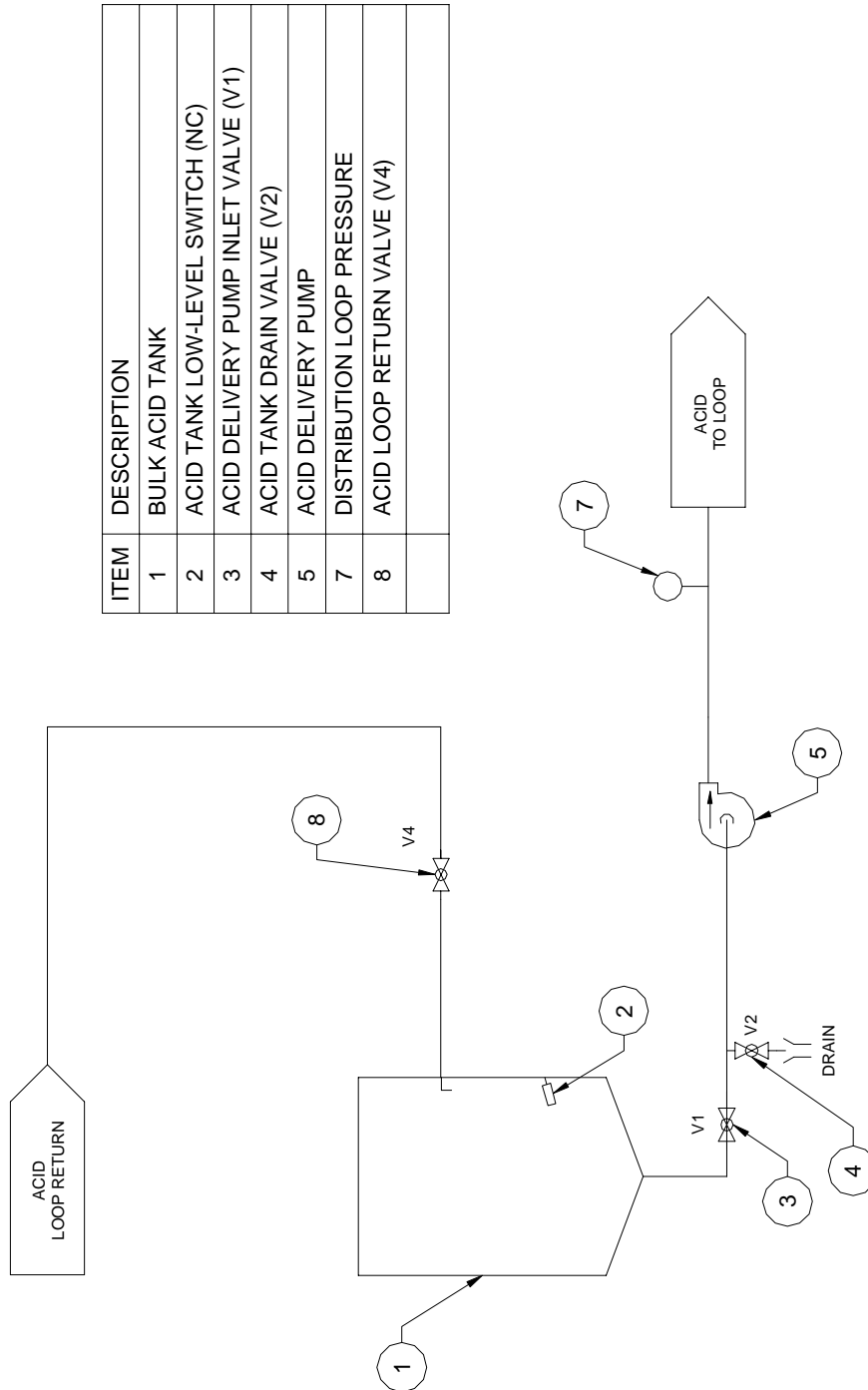

**Illustration 6-5: PSDS Controller Display Panel
(Three Acid Controller Shown)**

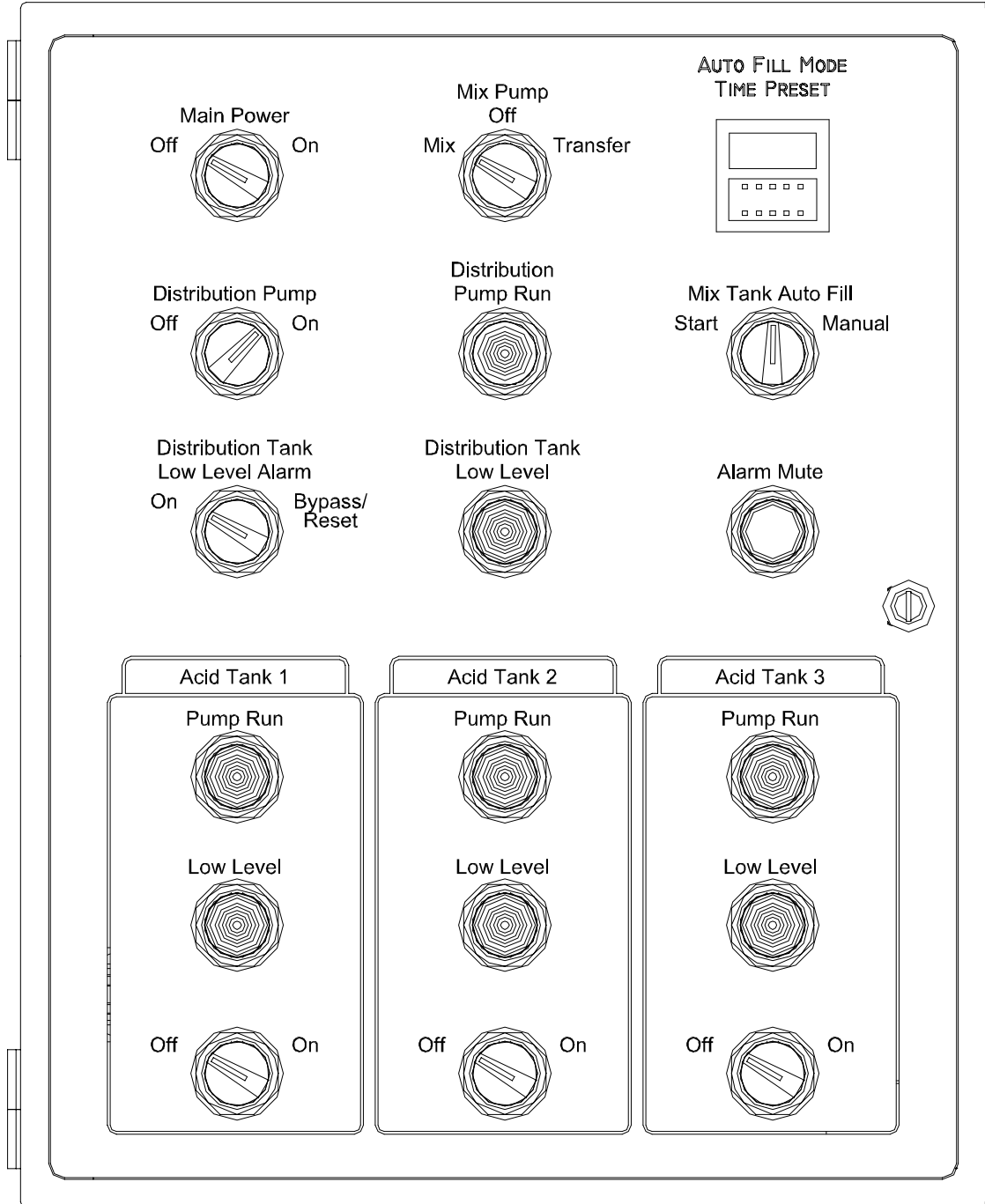

**Illustration 6-6: PSDS Controller Display Panel (Rear View)
(Three Acid Controller Shown)**




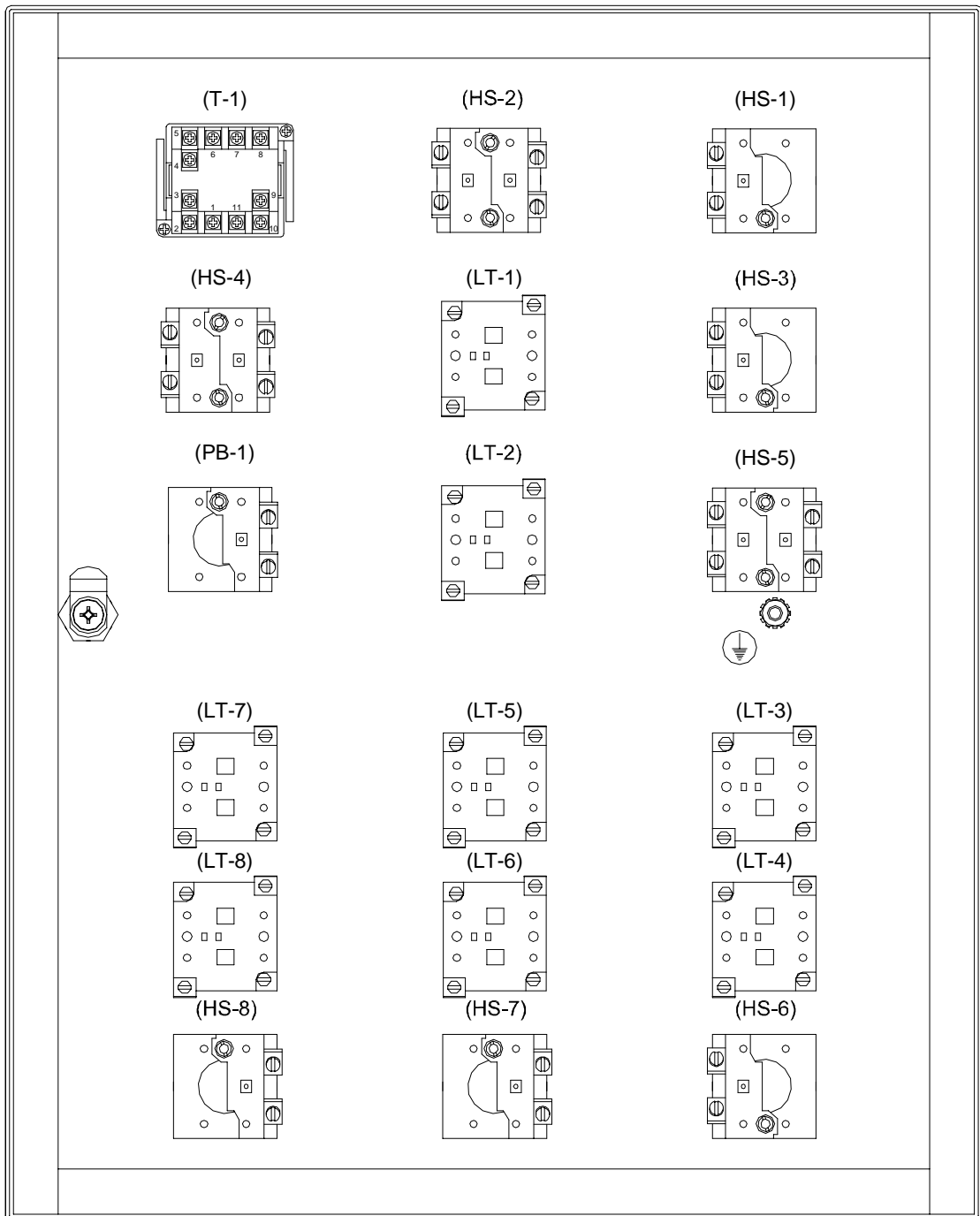


Illustration 6-8 VersaMax PLC

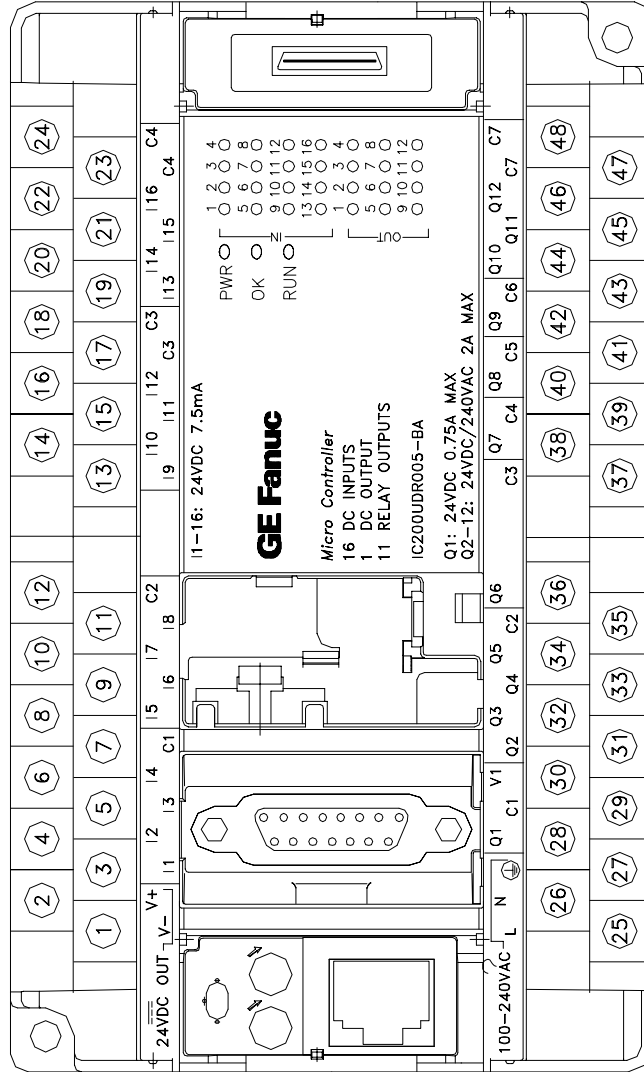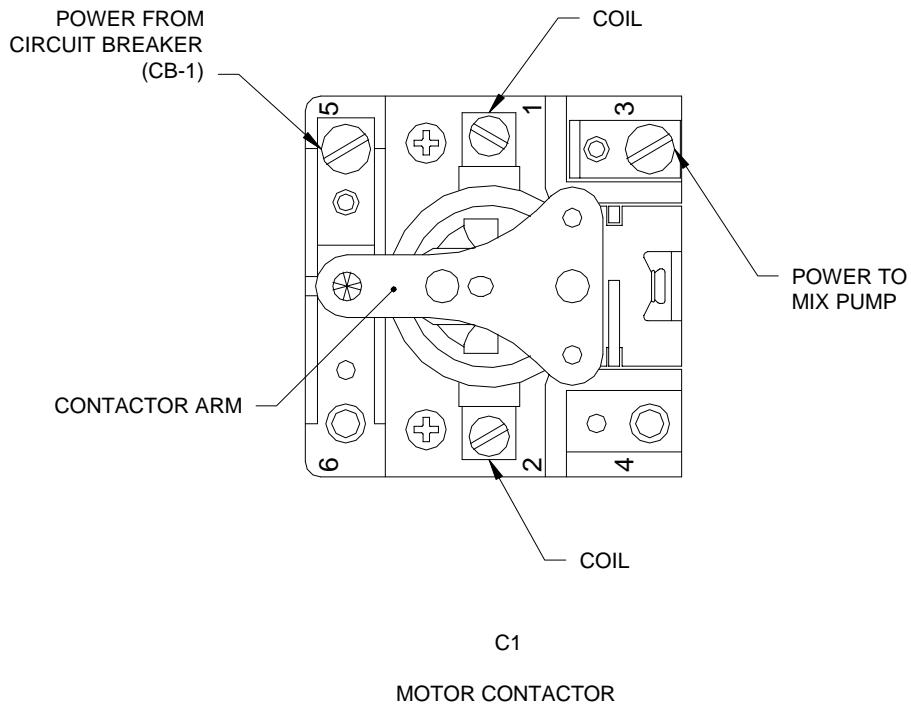

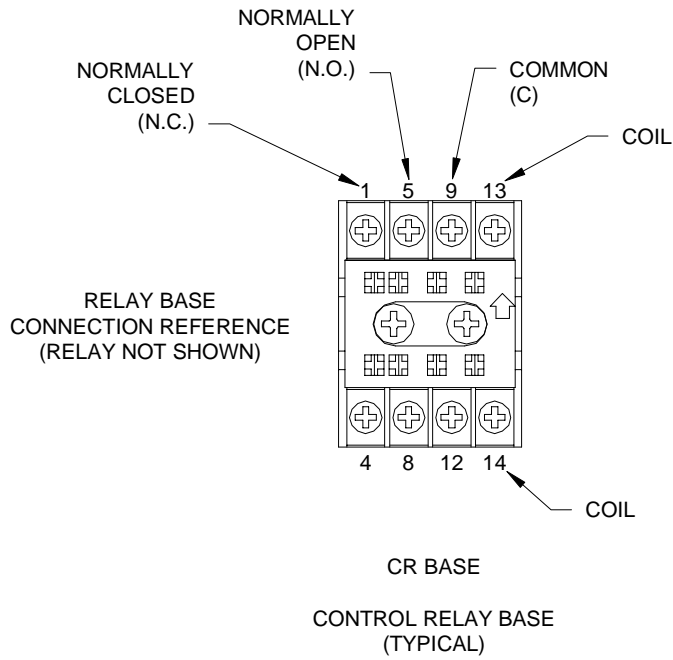

**Illustration 6-9: Electrical Components
(Control Relay and Motor Contactor)**

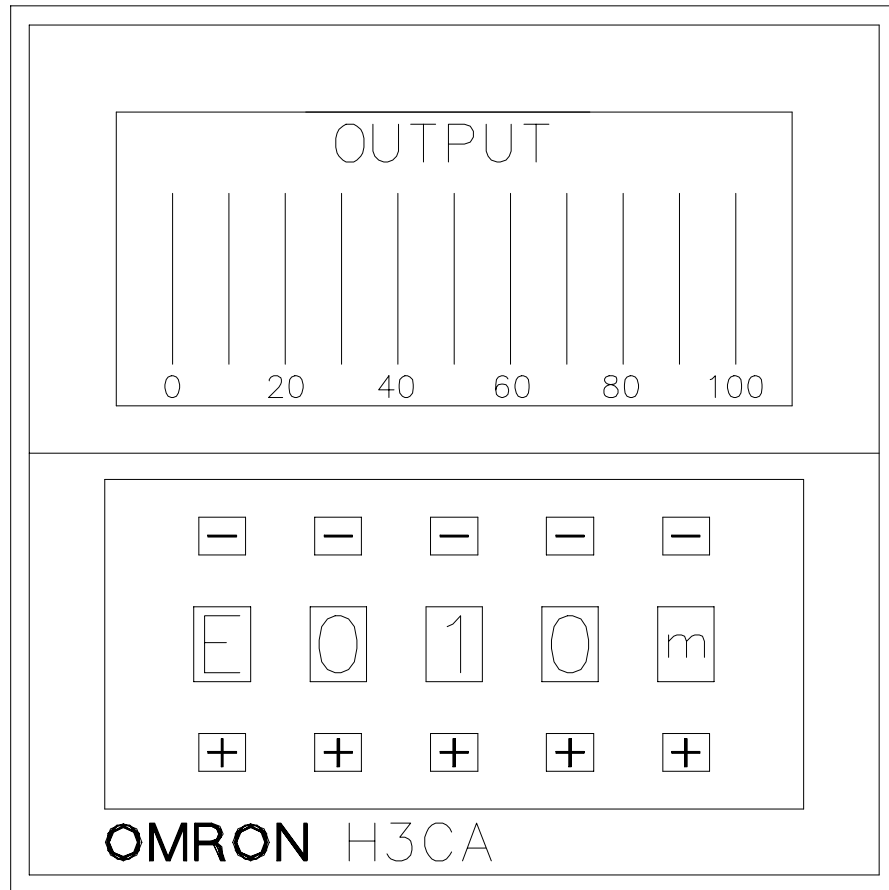


Illustration 6-10: Electrical Component (Auto Fill Timer)


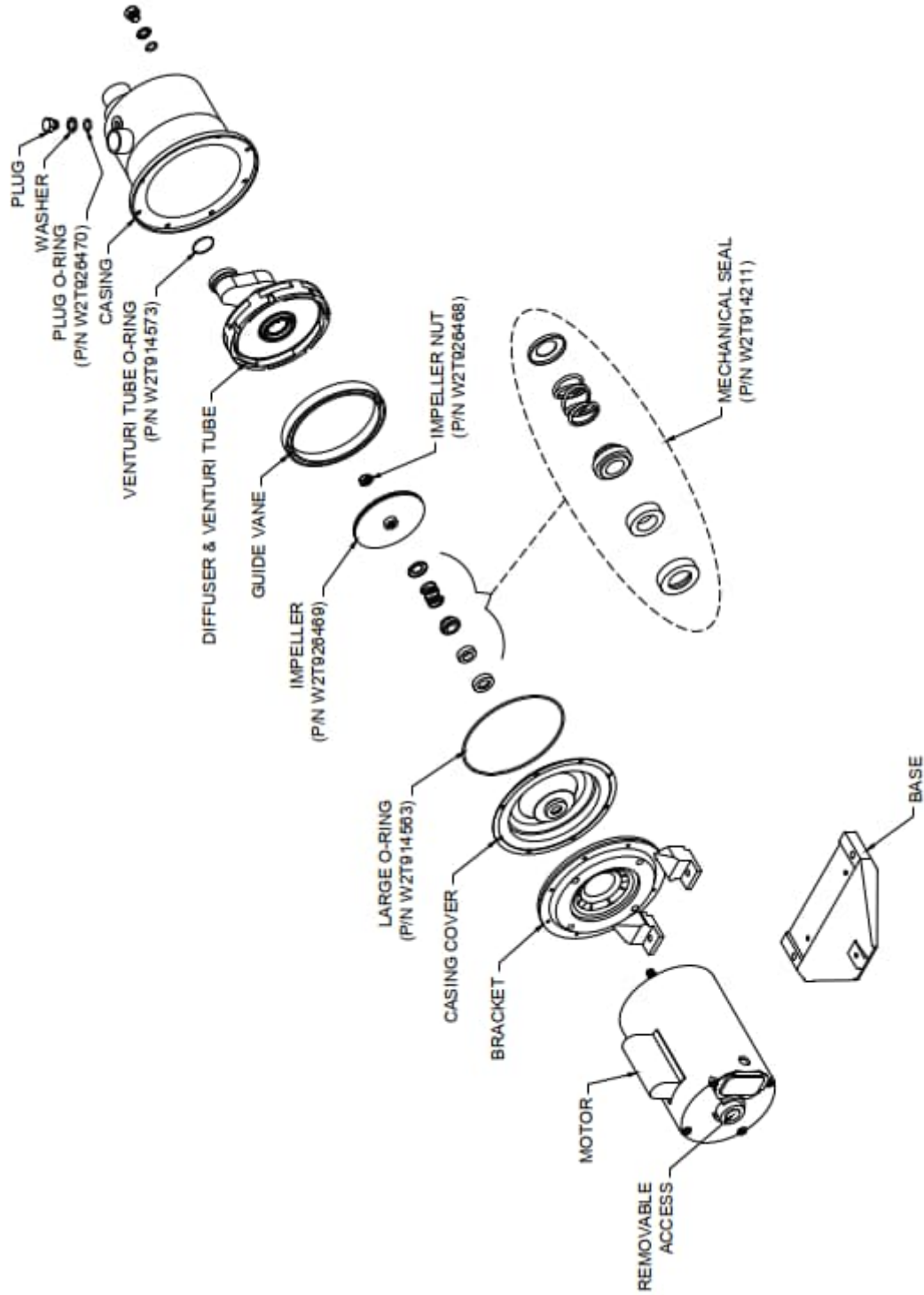


Illustration 6-11: Mix Pump/Motor Detail

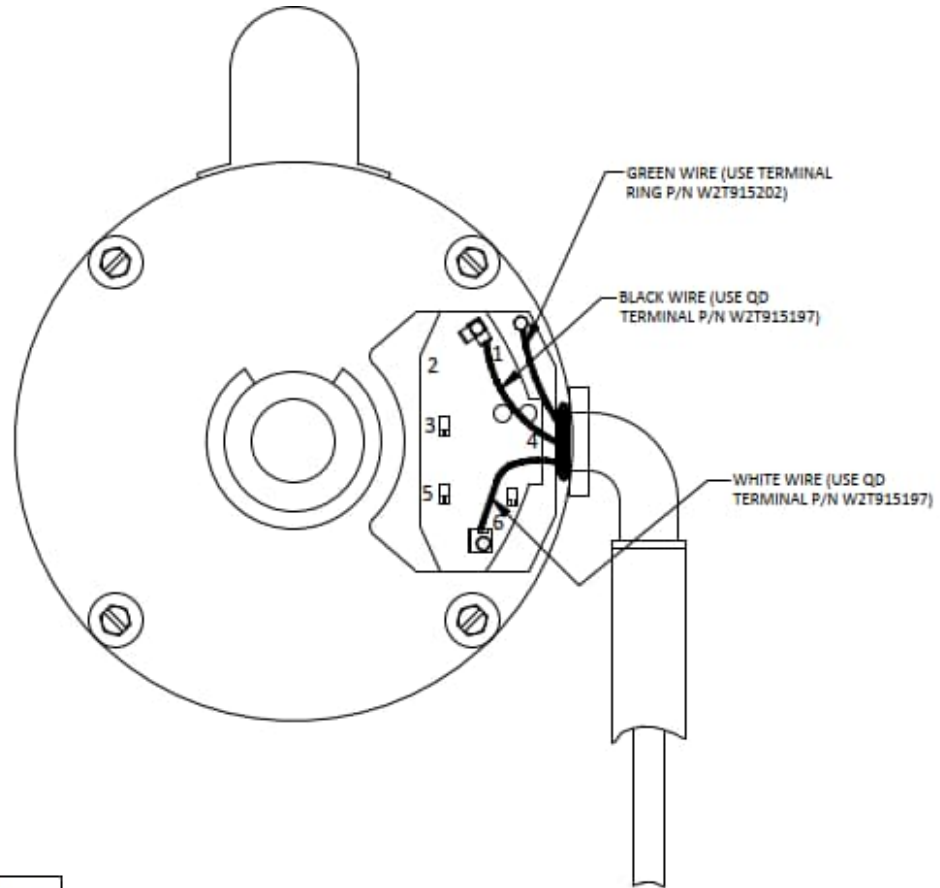

115 VAC WIRING

THE FOLLOWING ARE THE STEPS REQUIRED TO RECONFIGURE A REPLACEMENT MOTOR FROM 230 TO 115 VOLT OPERATION:

- 1) MOTOR CAN BE WIRED FOR 230 VOLT OR 120 VOLT OPERATION. THE SDS ELECTRICAL REQUIRES 115 VOLT CONFIGURATION.

WARNING: DISCONNECT ELECTRICAL POWER PRIOR TO SERVICING.

- 2) REMOVE COVER. DISCONNECT BROWN WIRE FROM TERMINAL MARKED 3 AND WHITE WIRE FROM TERMINAL MARKED 5.
- 3) RECONNECT WHITE WIRE TO TERMINAL MARKED 6.
- 4) RECONNECT BROWN WIRE TO TERMINAL MARKED 5.
- 5) CONNECT THE BLACK WIRE FROM 18-3 CORD (P/N W2T913166) TO TERMINAL 1. CONNECT THE WHITE WIRE TO TERMINAL 6, AND CONNECT THE GREEN WIRE TO THE GROUND SCREW.

WARNING: ROUTE WIRES AWAY FROM ROTATING CENTRIFUGAL SWITCH.

- 6) REPLACE COVER. REMOVE 230 VOLT STICKER FROM HOUSING. THIS WIRING IS FOR 115 VOLTS.

**Illustration 6-12: Mix Pump Electrical Configuration
(Rear of Motor Shown Without Cover)**

Illustration 6-13: PSDS Mix Tank

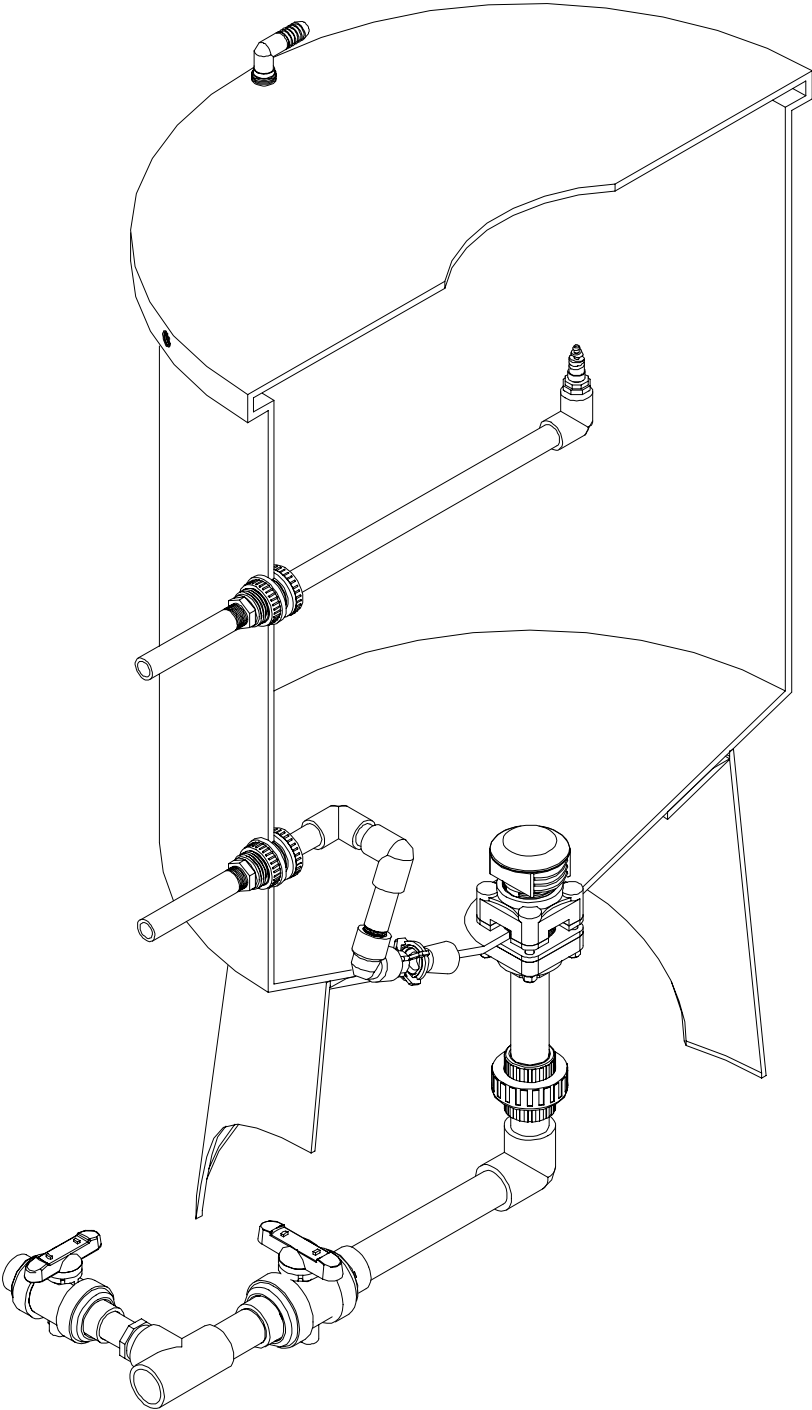


Illustration 6-14: PSDS Distribution Tank

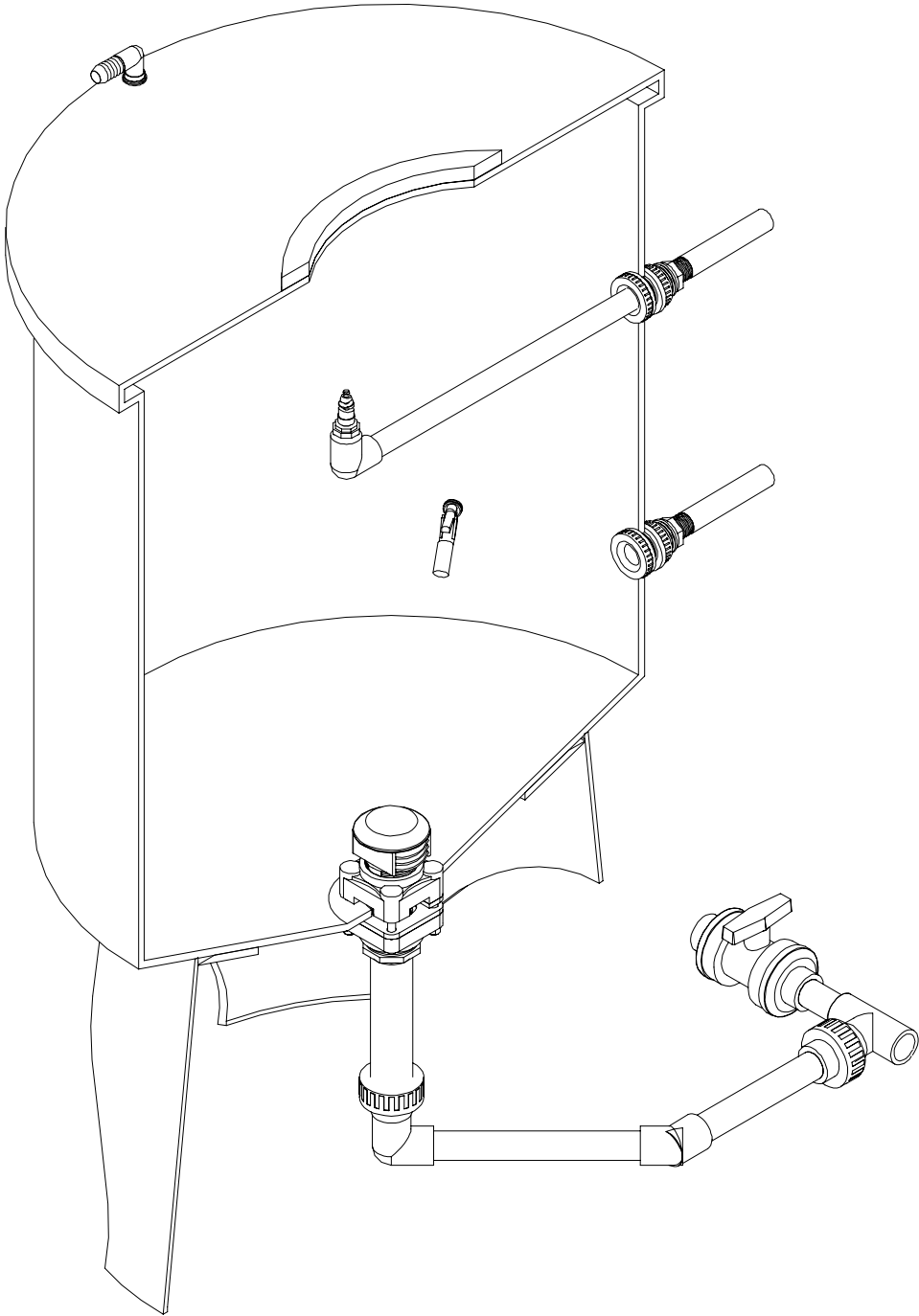


Illustration 6-15: PSDS Mix Tank To Frame Assembly

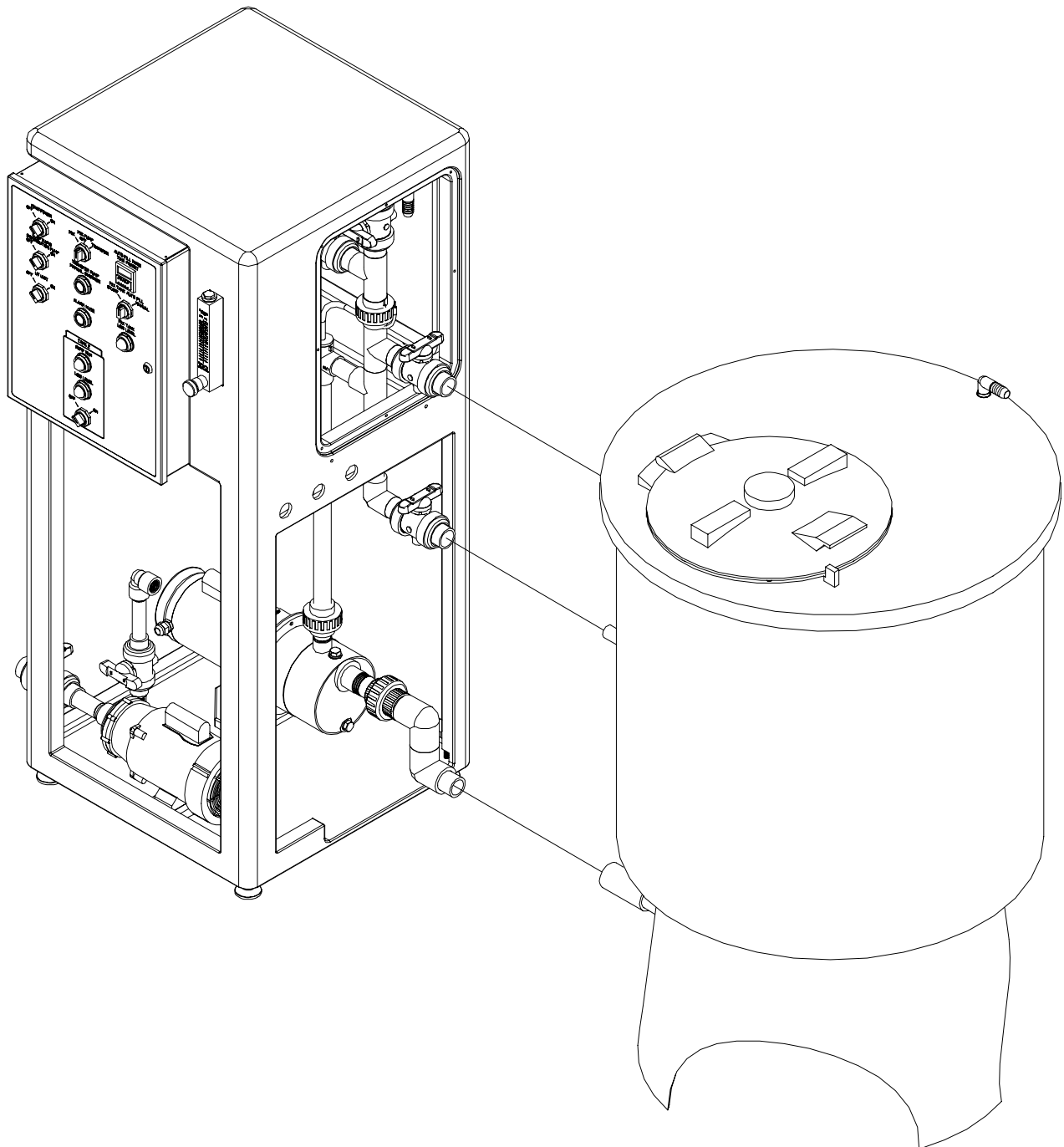


Illustration 6-16: PSDS Distribution Tank To Frame Assembly


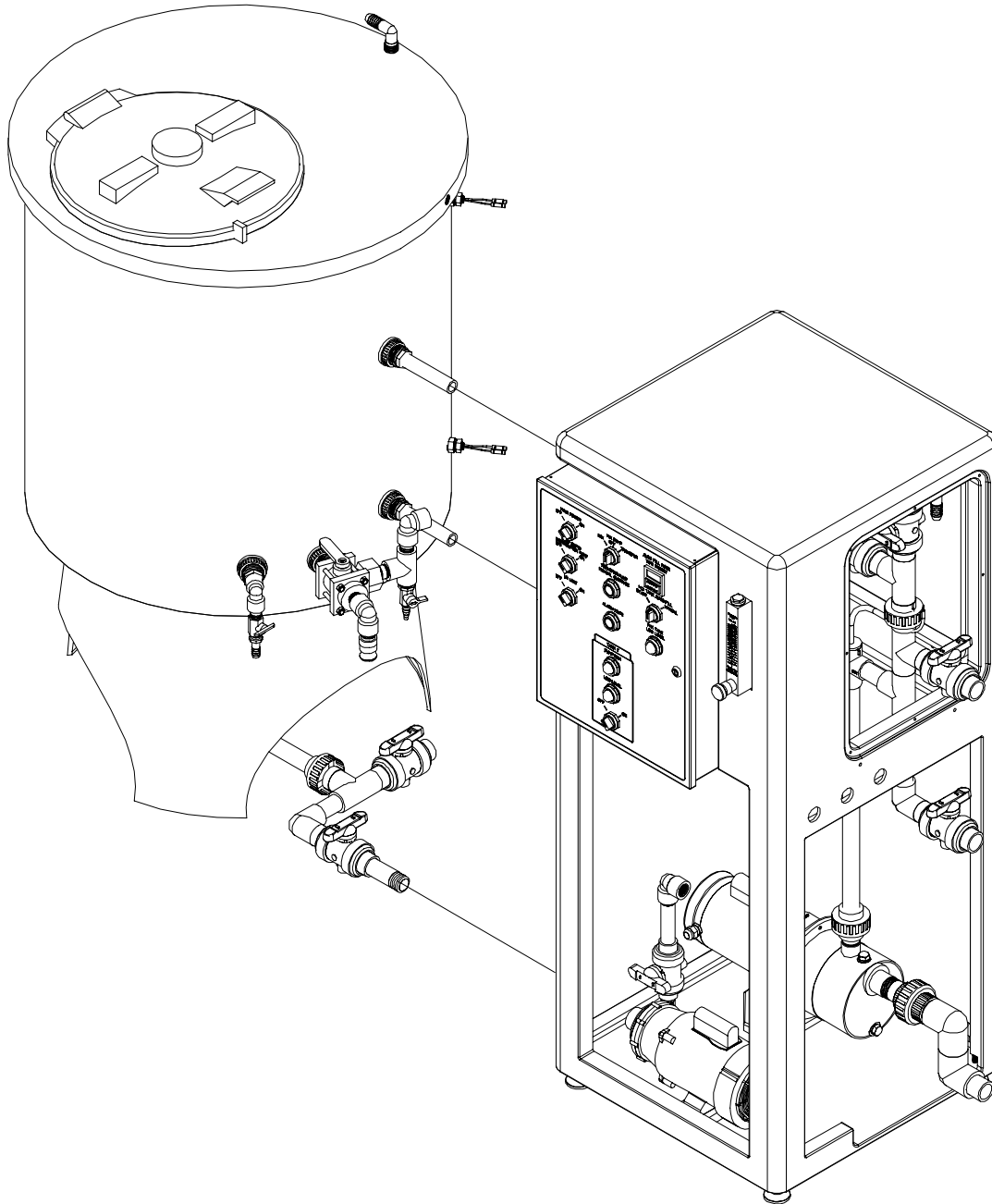


Illustration 6- 17: PSDS Drain Installation

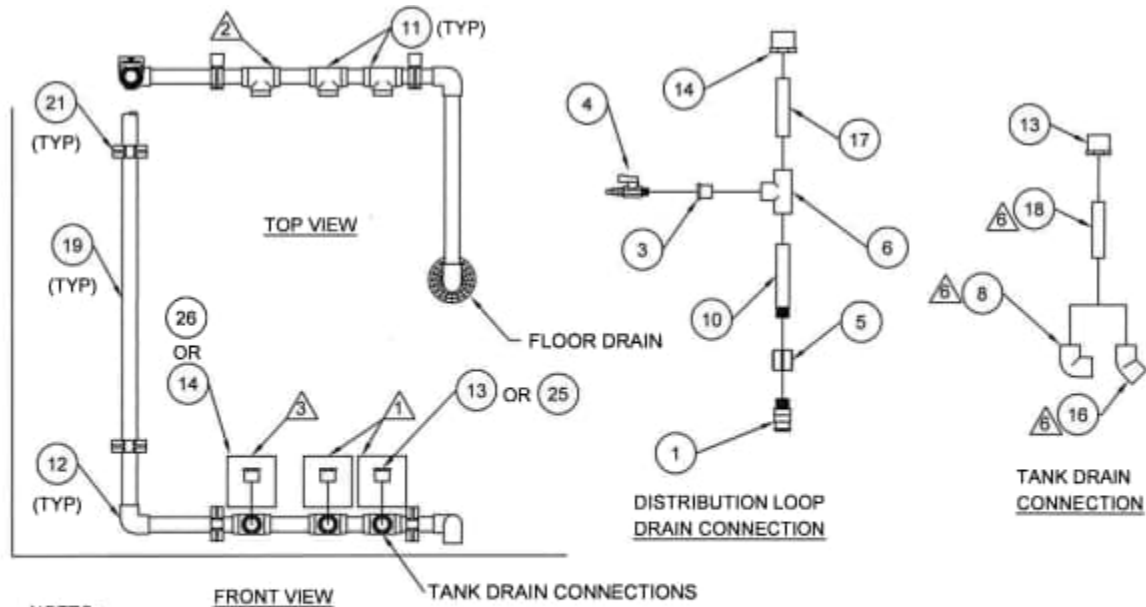


NOTES:

- ① INSERT THE 1-1/2" S X 1" THREADED BUSHING FOR THE TANK DRAIN CONNECTIONS.
- ② ANGLE TEE TO ALLOW EASE OF USE FOR GRAVITY FEED SAMPLE VALVE.
- ③ INSERT THE 1-1/2" S X 3/4" S REDUCER BUSHING FOR THE DISTRIBUTION LOOP DRAIN CONNECTION.
- ④ THE PSDS DRAIN CONNECTION ILLUSTRATIONS ARE SHOWN FOR REFERENCE PURPOSES EACH INSTALLATION IS DIFFERENT AND MAY REQUIRE VARIED CONNECTION CONFIGURATIONS.
- ⑤ DRILL AND TAP NPT HOLES AS REQUIRED FOR DRAIN CONNECTIONS. CONNECTIONS MAY BE PLACED IN THE VERTICAL OR HORIZONTAL DRAIN LINE.
- ⑥ USE WITH COMBINATION OF PARTS TO FIT AS NECESSARY.

26	BUSHING, RED, 1-1/2SX3/4T	W2T124957
25	BUSHING, RED, 1-1/2SX1T	W2T124958
24	CLAMP, WORM GEAR, 1", 5X445	W2T913584
23	QD, PP, 1"NPT, MALE	W2T913108
22	SCREW, DRYWALL, #6X1, PH, YELZINC	W2T919410
21	CLAMP, NYL, 1-1/2" PIPE, CLIC	W2T912719
20	WASHER, FLAT, SS, 1/4ID, APPROX. 65OD	W2T914026
19	PIPE, PVC, 1-1/2", SCH80	W2T913164
18	PIPE, PVC, 1", SCH80	W2T913162
17	PIPE, PVC, 3/4", SCH80	W2T913161
16	ELBOW, PVC, 1SXS, 45DEG, SCH80	W2T915695
15	ELBOW, 3/4"SLIP, 45DEG, SCH80	W2T914325
14	BUSHING, RED, 1-1/2SX3/4S	W2T219673
13	BUSHING, RED, 1-1/2SX1S	W2T914308
12	ELBOW, 1-1/2SXS, SCH80	W2T206239
11	TEE, PVC, 1-1/2" SX 1-1/2" SX 1-1/2" S, SCH 80	W2T356688
10	NIPPLE, PVC, 3/4X12, SCH80	W2T914501
9	ELBOW, PVC, 3/4", SXS	W2T219471
8	ELBOW, PVC, 1SXS, SCH80	W2T494919
7	COUPLER, PVC, 1-1/2SXS, SC80	W2T914251
6	TEE, PVC, 3/4SXSXS, SCH 80	W2T914301
5	COUPLER, PVC, 3/4XTX, SCH80	W2T913091
4	VALVE, BALL, 1/4MNPTXBARB	W2T913086
3	BUSHING, RED, 3/4SX1/4T	W2T124949
2	COUPLER, PVC, 3/4SXS, SCH80	W2T127324
1	QD, PP, 3/4"NPT, MALE	W2T914862
ITEM	DESCRIPTION	PART NO.


Illustration 6-18: Bicarb Distribution Loop Connections

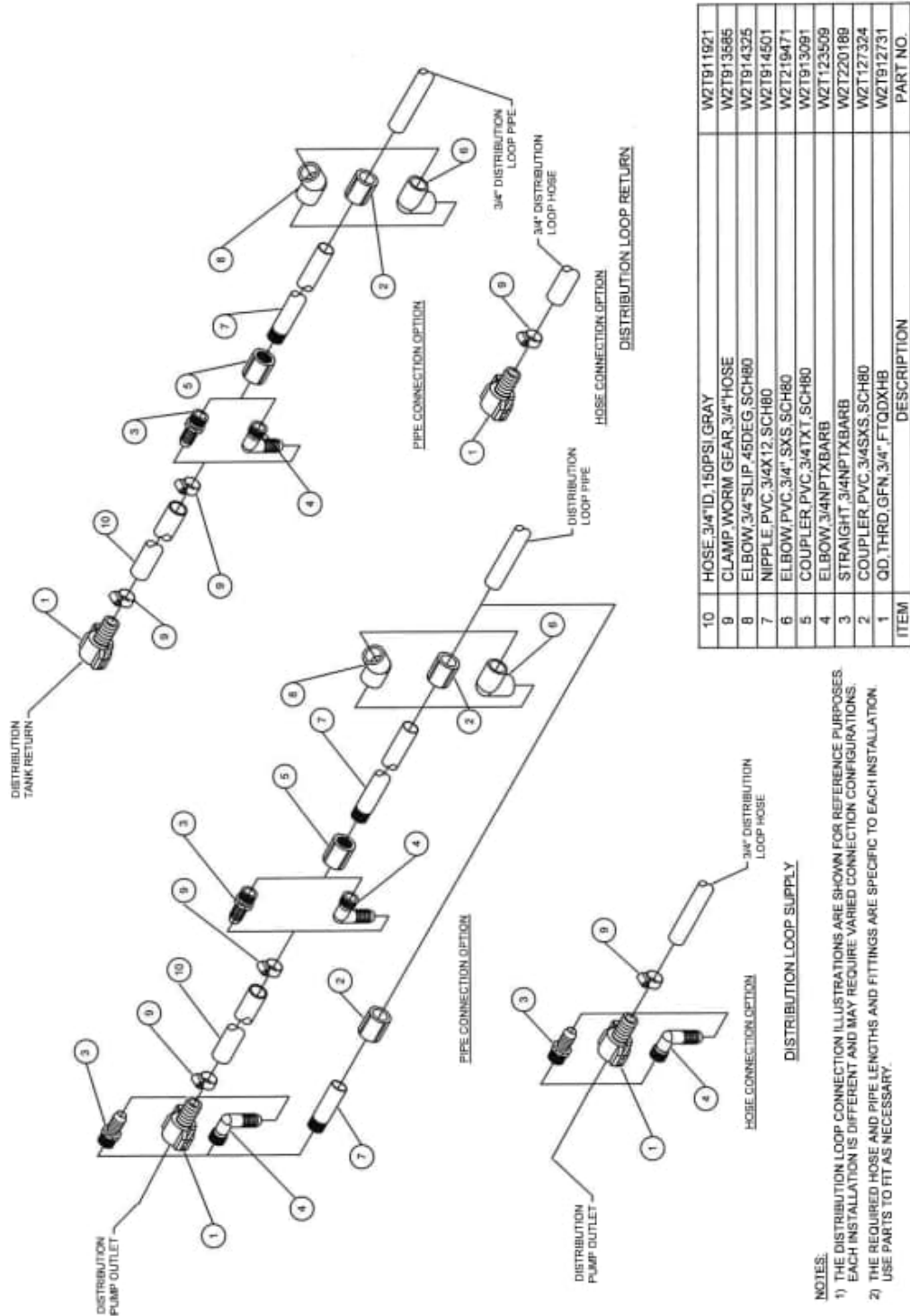

NOTES:
 1) THE DISTRIBUTION LOOP CONNECTION ILLUSTRATIONS ARE SHOWN FOR REFERENCE PURPOSES. EACH INSTALLATION IS DIFFERENT AND MAY REQUIRE VARIOUS CONNECTION CONFIGURATIONS.
 2) THE REQUIRED HOSE AND PIPE LENGTHS AND FITTINGS ARE SPECIFIC TO EACH INSTALLATION. USE PARTS TO FIT AS NECESSARY.

ITEM	DESCRIPTION	PART NO.
10	HOSE, 3/4" ID, 150PSI, GRAY	W2T911921
9	CLAMP, WORM GEAR, 3/4" HOSE	W2T913585
8	ELBOW, 3/4" SLIP, 45DEG, SCH80	W2T914325
7	NIPPLE, PVC, 3/4X12, SCH80	W2T914501
6	ELBOW, PVC, 3/4", SXS, SCH80	W2T219471
5	COUPLER, PVC, 3/4" TX, SCH80	W2T913091
4	ELBOW, 3/4" NPT X BARB	W2T123509
3	STRAIGHT, 3/4" NPT X BARB	W2T220189
2	COUPLER, PVC, 3/4" SXS, SCH80	W2T127324
1	OD, THRD, GFN, 3/4", FTQDXHB	W2T912731

Illustration 6-19: Bicarb Distribution Loop Pump Connections

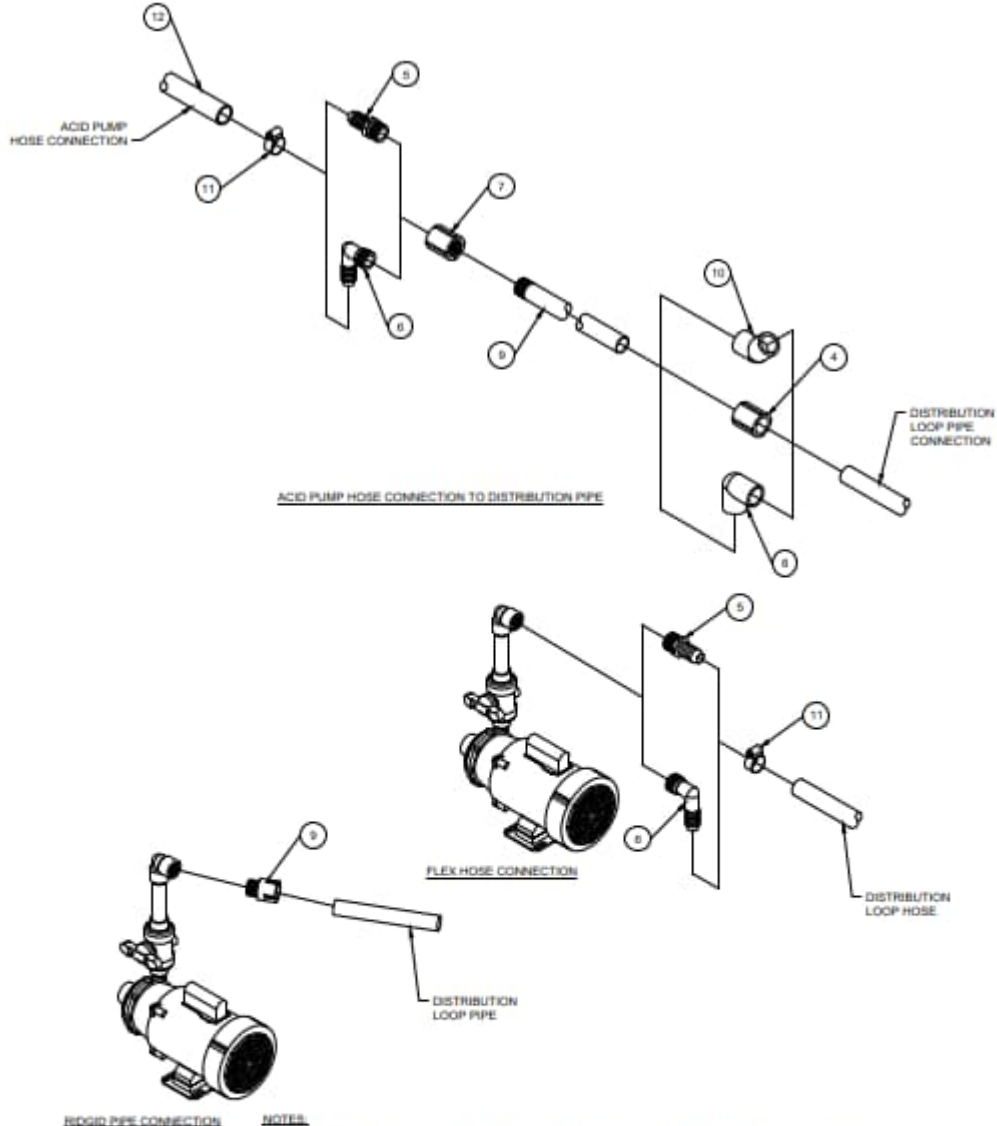

RIGID PIPE CONNECTION

NOTES:

- 1) THE DISTRIBUTION LOOP CONNECTION ILLUSTRATIONS ARE SHOWN FOR REFERENCE PURPOSES. EACH INSTALLATION IS DIFFERENT AND MAY REQUIRE VARIED CONNECTION CONFIGURATIONS.
- 2) THE REQUIRED HOSE AND PIPE LENGTHS AND FITTINGS ARE SPECIFIC TO EACH INSTALLATION. USE PARTS TO FIT AS NECESSARY.

12	HOSE, 3/4"ID, 150PSI, GRAY	W2T911921
11	CLAMP, WORM GEAR, 3/4"HOSE	W2T913585
10	ELBOW, 3/4SLIP, 45DEG, SCH80	W2T914325
9	NIPPLE, PVC, 3/4X12, SCH80	W2T914501
8	ELBOW, PVC, 3/4"XKS, SCH80	W2T219471
7	COUPLER, PVC, 3/4TXT, SCH80	W2T913091
6	ELBOW, 3/4NPTXBARB	W2T123509
5	STRAIGHT, 3/4NPTXBARB	W2T220189
4	COUPLER, PVC, 3/4SXS, SCH80	W2T127324
3	ADAPTER, PVC, 3/4SXT	W2T124737
2	PUMP DISTRIBUTION, PSDS	W2T914732
1	PUMP, MAG, PP, 1/3HP, 115/230, TE-5 SC-MD, MAR	W2T911979
ITEM	DESCRIPTION	PART NO.

Illustration 6-20: PSDS Acid Tank Connections

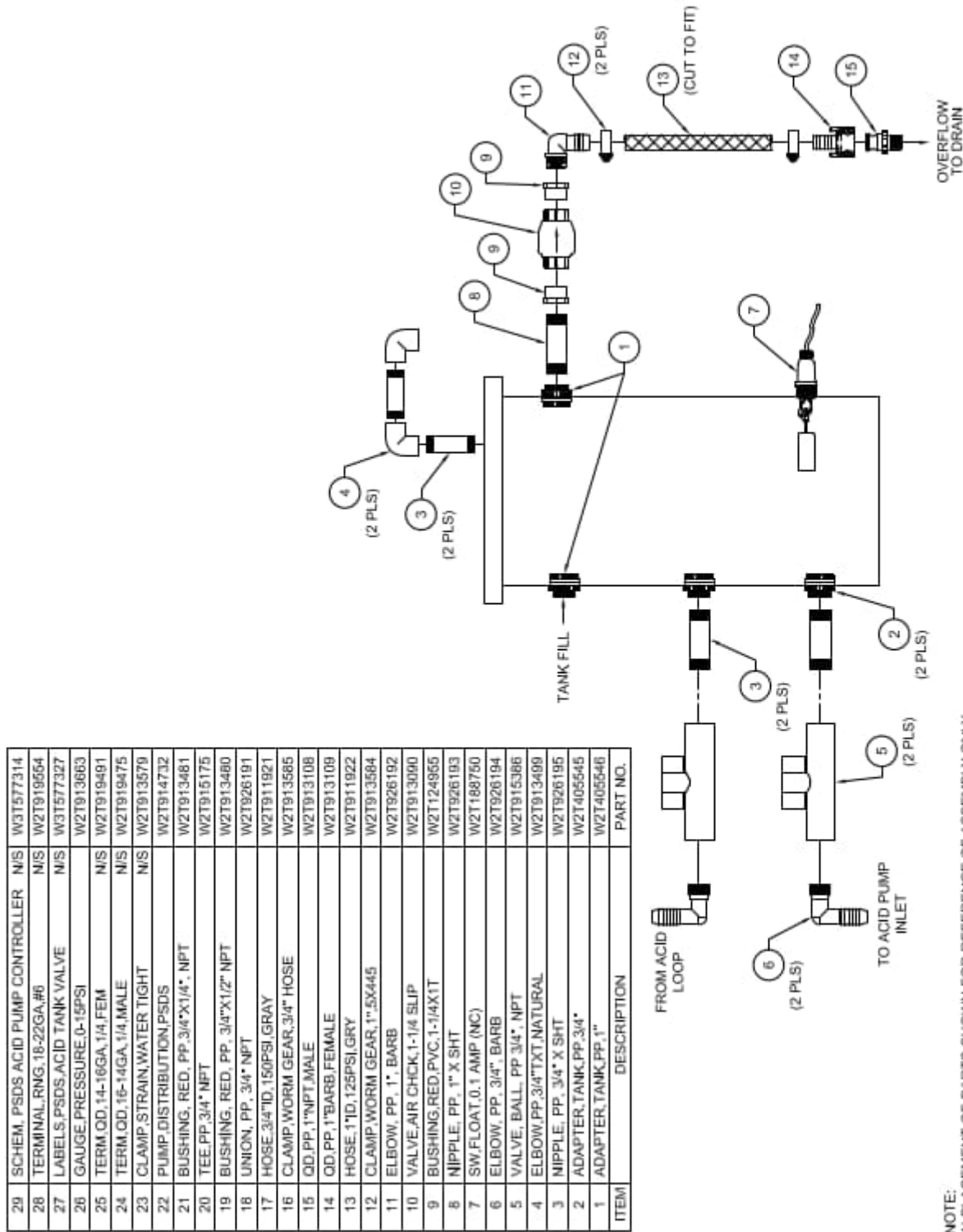


Illustration 6-21: PSDS Status Monitor, Surface Mount Models



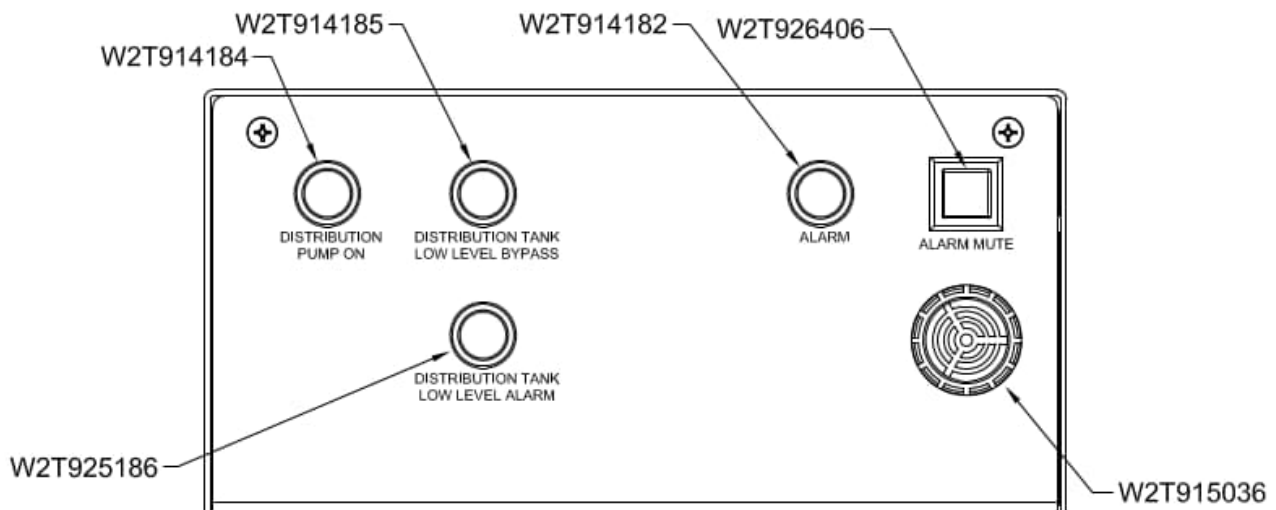


Illustration 6-22: PSDS Status Monitor, Flush Mount Models

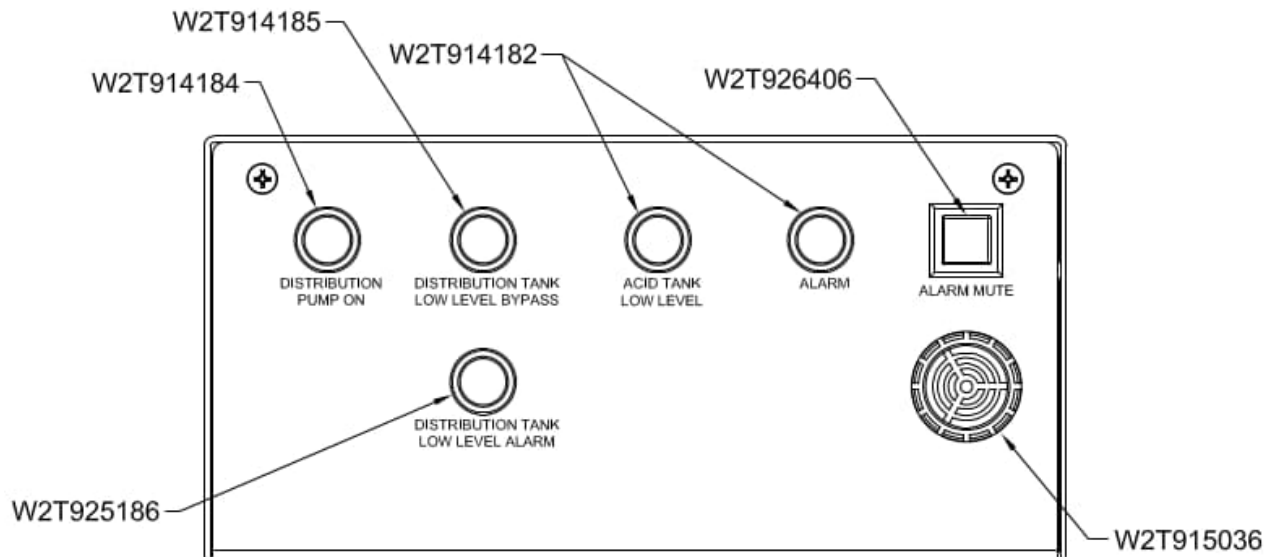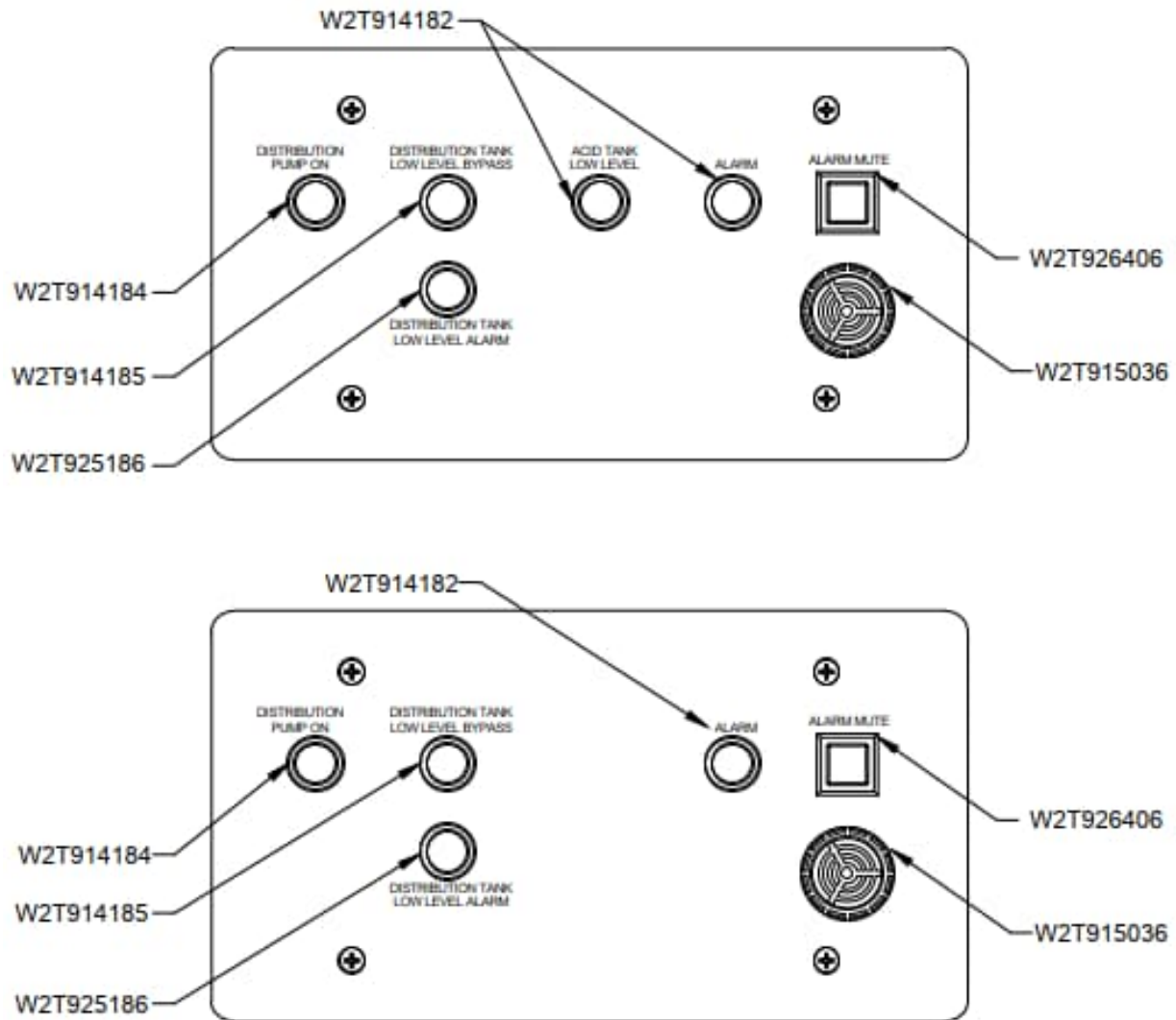

W2T925186	INDICATOR, LED, 1/4", 24VDC, RED
W2T926406	SWITCH, PUSH BUTTON, MOM ON, BLACK
W2T914185	INDICATOR, LED, 1/2", 24V, AMBER
W2T915036	SONALERT, 4-28VDC, PNL MNT
W2T914184	INDICATOR, LED, 1/2", 24V, GREEN
W2T914182	INDICATOR, LED, 1/2", 24V, RED
ITEM	DESCRIPTION

NOTES:

**Pressurized Solution Distribution System
PSDS**

CHAPTER SEVEN:
SPARE PARTS LIST

This page intentionally left blank.

PSDS SPARE PARTS LIST

CATALOG NUMBER	DESCRIPTION
W2T913654	Fuse, 1 Amp Slo Blo
W2T913070	Fuse, 1.5 Amp Slo Blo
W2T913457	Fuse, 5 Amp, Slo Blo
W2T913866	Filter, Carbon, 9.75", 5700
W3T578377	Wire, Dip Tube, 45'
W3T579417	PLC, Versa Max Main Module
W2T188750	Switch, Float
W2T914910	Relay, 24 VAC
W2T914915	Relay, 24 VAC, 20 Amp (MCR-1)
W2T915035	Sonalert, 24 Volt
W2T915248	Timer, Digital
W2T914918	Relay, 24 VAC, 30 Amp (C-1)
W2T919762	Lamp, Red, 24 Volt
W2T919761	Lamp, Green, 24 Volt
W2T913554	Circuit Breaker, 20 Amp
W2T914862	QD, Male, 3/4" MNPT, Cam-Loc
W2T448912	QD, Female, 3/4" Barb, Cam-Loc
W2T393358	Valve, Ball, Union, 1"
W2T913086	Valve, Ball, 1/4" X Barb
W2T369988	Valve, Ball, Union, 3/4"
W2T892384	Valve, Ball, Union, 1/2"
W2T915059	Spray Head, PVC, 3/8", Mix Tank
W2T913108	QD, Male, 1" MNPT, Cam-Loc
W2T913109	QD, Female, 1" Barb, Cam-Loc
W2T913116	Switch, Flow, 1", 2 GPM
W2T913122	Valve, Solenoid, Plastic (Water Inlet)
W2T362196	Flow Meter, W/Needle, 5 GPM
W2T913763	Eductor, Mixing
W2T919610	Tube, Coil, PVC, Clear
W2T913156	Pump, 3/4 HP, Mix Pump, EBARA
W2T914997	Silicone Sealant, 3 oz
W2T926470	O-ring, Vent Plug Ebara
W2T914573	O-ring, Venturi Tube
W2T914563	O-ring, Pump Case, Ebara
W2T914211	Kit, Seal, Shaft, Ebara
W2T926469	Impeller, Ebara Pump
W2T926468	Nut, Ebara Impeller
W2T924570	Tank Lid (Hinged), Non-Vented
W2T425370	Switch, 3 Pos, Maintain, 2NO, 2NC, GE

CATALOG NUMBER	DESCRIPTION
W2T38399	Switch, P-Button, 1NO, GE
W2T926415	Switch, 3 Pos, Mom, 1NO, 1NC, GE
W2T425373	Switch, 2 Pos, Maintain, 1NO, 1NC, GE
W2T926413	Indicator, 24V Red, GE
W2T911932	Indicator, 24V Green, GE
W2T913755	Drain, PSDS
W2T913754	Drain, Anti-Vortex
W2T915108	Switch, 2 Pos, 24v, Illuminated, Amber
W2T800311	Indicator, 120V, Green
W2T911968	Fuse, 2.5 Amp Slow Blow
W2T914554	Hinged Lid O-ring
W2T913022	Valve, Ball, Union, 1-1/4"
W2T122525	Valve, Ball, 3-Way, 3/4" Threaded
W2T911979	Pump, Distribution
W2T913663	Pressure Gauge, 15 PSI
W2T913083	Valve, Ball, 1/4" FNPT x MNPT

PSDS STATUS MONITOR SPARE PARTS LIST

CATALOG NUMBER	DESCRIPTION
W2T914182	Indicator, LED, 1/2", 24VAC, Red
W2T914184	Indicator, LED, 1/2", 24VAC, Green
W2T914185	Indicator, LED, 1/2", 24VAC, Yellow
W2T925186	Indicator, LED, 1/2", 24VDC Red
W2T915036	Sonalert, 4-28VDC, Pnl Mnt
W2T913579	Clamp Strain, Water Tight
W2T926406	Pushbutton Switch
W2T915193	Terminal, Butt Splice, No-Crimp

**Pressurized Solution Delivery System
PSDS**

**APPENDIX A:
TECHNOTES/LOG SHEETS**

This page intentionally left blank.

TechNotes

In this section, you will find a variety of information related to your equipment and to dialysis water use in general.

Extra copies of TechNotes may be ordered from Customer Service at (800) 633-3080.

Please reference the number of the Technote, preceded by the letters TN; i.e., Technote 133 = TN133.

This page intentionally left blank.

Call **(800-633-3080)** for additional information or visit **www.mcpur.com**.



14550 28th AVE N
Plymouth, MN
55447
USA
Ph: 800-633-3080
FAX: 763-210-3868