

Tx3[®]

THERAPY SUPPLY STATION



MAR COR

PURITY FOR LIFE

A reimagined dialysis workstation that's clean and durable. As connection stations do not separate product water ports and acid & bicarbonate concentrates from waste solution ports, the risk of microorganism exposure can be increased. All aspects of the Tx3 Therapy Supply Station are focused on maintaining a clean and safe environment. Unlike other wall boxes, the Tx3 station is:

VERSATILE Accommodates a variety of water and port options, including high purity and standard water connections.

ADAPTABLE Fits current model cut-outs for easy retrofits.

COMPLIANT Meets the CDC recommended strategy of separate areas, with the introduction of a door compartment that separates the waste drain, minimizing the risk of cross-contamination.

IDENTIFIABLE Reduces staff confusion with clear labels and color-coded ports and valve handles, meeting AAMI and CMS standards.

CLEAN Easy to keep clean with smooth surfaces, a removable face cover and sloped areas that also aid in complete fluid drainage.

DURABLE Made of corrosion resistant materials for a dependable, long life.

Tx3[®] THERAPY SUPPLY STATION

Materials of Construction:

ABS resin

Water Connection:

PEX, PFA, PVDF, and other standard connections

Concentrate Connector:

Fresenius or B.Braun

Drain Connection:

1-1/2" male socket

Mounting Kit:

Included with each station (Stainless steel screws and washers)

Dimensions:

Height: 17.9" (19.3" with bezel)

Width: 17.3" (18.2" with bezel)

Depth: 5" (6" with bezel)

Weight: 3.5 lbs. (empty) 4.5 lbs. (with bezel)

Rough-in: 15" W x 15-3/4" H

Mounting face area: 1" minimum




Configuration Examples

CleanPEX™ Tubing Loops for Water

- **DESCRIPTION:** Wall Box, Tx3, WP, FP, 2A, 1B, 25MM
- **DETAIL:** CWP water fitting (25MM PEX). Fresenius connectors, 2 acid and 1 bicarb (color coded). Shut-off valve per connector.

PFA Tubing Loops for Water

- **DESCRIPTION:** Wall Box, Tx3, WT, FP, 2A, 1B, 0.75
- **DETAILS:** Water fitting (3/4" PFA). Fresenius connectors, 2 acid and 1 bicarb (color coded). Shut-off valve per connector.

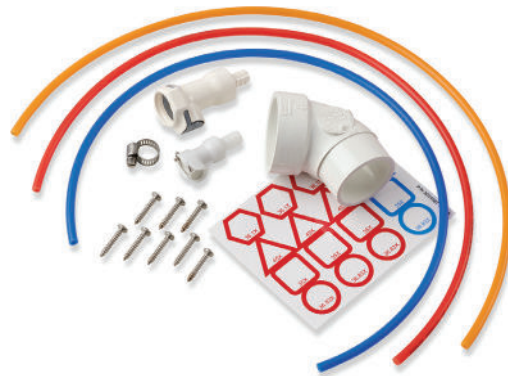

Visit mcpur.com for more information or call 1-800-633-3080

Tx3 is a trademark of Evoqua Water Technologies LLC, its subsidiaries or affiliates in some countries. Tx3 Therapy Supply Station is patented in some countries. CleanPEX is a trademark of Uponor Innovation AB.

All information presented herein is believed reliable and in accordance with accepted engineering practices. Evoqua makes no warranties as to the completeness of this information. Users are responsible for evaluating individual product suitability for specific applications. Evoqua assumes no liability whatsoever for any special, indirect or consequential damages arising from the sale, resale or misuse of its products.